
# User manual

Addax

**R82**<sup>®</sup>  
by Etac


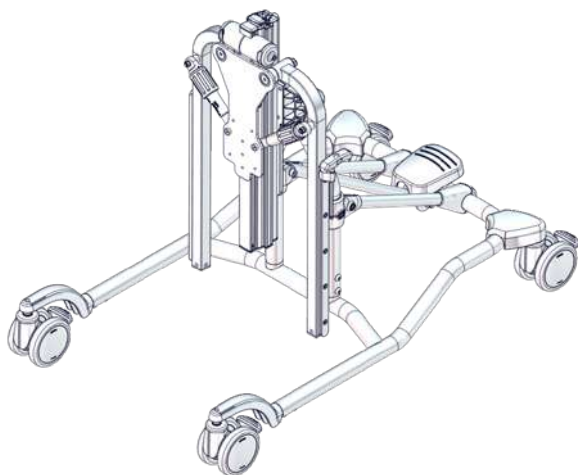
M1622 - 2025-12-16 - Rev. 5.0

[www.etac.com](http://www.etac.com)

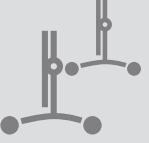
Content

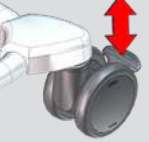	Technical information	9
		10-15
		16
mm (inch) kg (lb)		17-19
		20
		21-25
		26


---

## Content


---

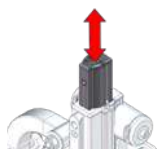
27

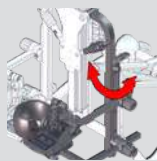28

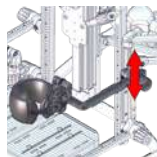29

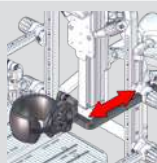
29

30


30

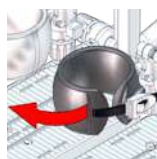
31

---


## Content

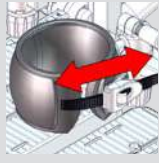
31

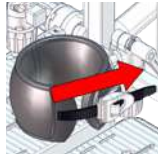
32

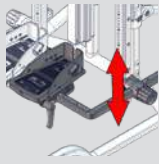
33

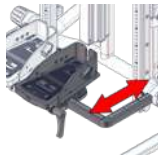
33

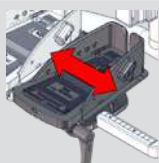34

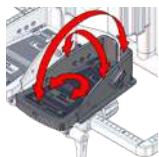34

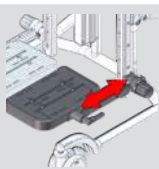
35


---

## Content

---

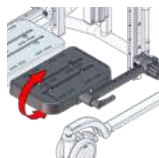
35

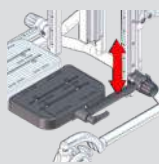36


37

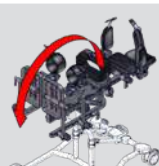38

39

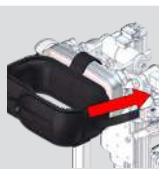40


40

---


## Content

41

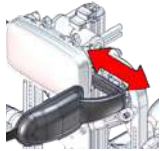
42

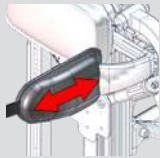42-43


44

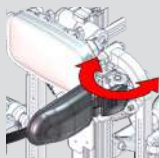44-45


45

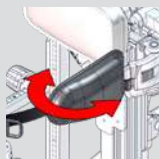
46

# Content

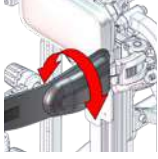46

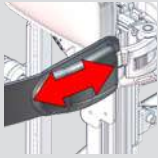47

48-53

54-55

56

57-59

---

## LANGUAGE

ENGLISH (Original instructions)	60-63
DEUTSCH	64-69
NEDERLANDS	70-75
DANSK	76-79
NORSK	80-83
SVENSKA	84-87
SUOMI	88-91
ESPAÑOL	92-96
FRANÇAIS	97-102
ITALIANO	103-107
PORTUGUÊS	108-113
PORTUGUÊS DO BRASIL	114-118
РОССИЯ	119-124
POLSKI	125-130
ČESKY	131-134
中文	135-138
日本語	139-144
ΕΛΛΗΝΙΚΑ	145-150
ÍSLENSKA	151-154

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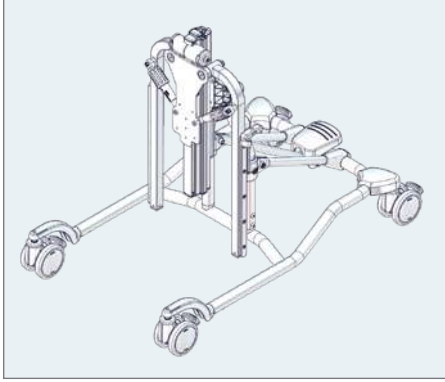
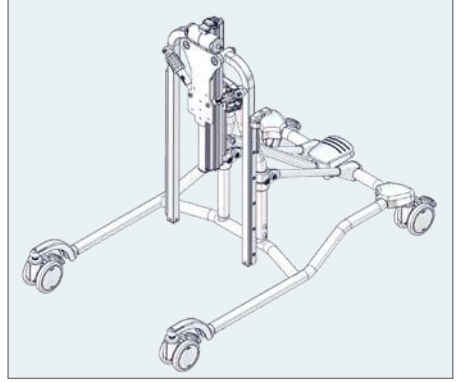
**Size 2**

## Size 1 & 2

		mm (inch), kg (lb)	
		Size 1 	Size 2 
	Maximum user mass	30 kg (66 lb)	45 kg (99 lb)
	User height	700 - 1050 mm (27½ - 41¼")	930 - 1280 mm (36½ - 50½")
	Foot supports, parted to axil height	370 - 800 mm (14½ - 31½")	370 - 1060 mm (14½ - 41¾")
	Sandals, 3D adjustable, to axil height	330 - 760 mm (13 - 30")	330 - 1020 mm (13 - 40¾")
	Leg length from foot supports, parted to pivot point	170 - 500 mm (6¾ - 19¾")	170 - 650 mm (6¾ - 25½")
	Leg length from sandals, 3D adjustable to pivot point	130 - 460 mm (5 - 18")	130 - 610 mm (5 - 24")
	Tray height from foot supports, parted, in supine	170 - 720 mm (6¾ - 28¾")	170 - 980 mm (6¾ - 38½")
	Tray height from sandals, 3D adjustable, in supine	130 - 760 mm (5 - 30")	130 - 940 mm (5 - 37")
	Tray height from foot supports, parted, in prone	290 - 800 mm (11½ - 31½")	290 - 1070 mm (11½ - 42¾")






mm  
(inch)  
kg  
(lb)







Size 1 & 2			
		mm (inch), kg (lb)	
		Size 1 	Size 2 
	Tray height from sandals, 3D adjustable, in prone	250 - 760 mm (9¾ - 30")	250 - 1030 mm (9¾ - 40½")
	Entrance height, lowest height, for foot supports, parted	145 mm (5¾")	145 mm (5¾")
	Entrance height, lowest height, for sandals, 3D adjustable	185 mm (7¼")	185 mm (7¼")
	Abduction angle, each leg	0 - 30 °	0 - 30 °
	Foot supports, parted angle	(-30) - 30 °	(-30) - 30 °
	3D adjustment angle, sandals	25 °	25 °
	Tilt	0 - 90 °	0 - 90 °
	Support surface height, horizontal	760 mm (30")	840 mm (33")
	Overall height	1100 - 1480 mm (43¾ - 58¼")	1230 - 1610 mm (48½ - 63½")

mm  
(inch)  
kg  
(lb)

## Size 1 & 2


		mm (inch), kg (lb)	
		Size 1 	Size 2 
	Overall length	1060 - 1490 mm (41¾ - 58¾")	1200 - 1620 mm (47¼ - 63¾")
	Overall width	540 - 1030 mm (21¼ - 40½")	620 - 1140 mm (24½ - 45")
	Static stability	10 °	10 °
	Mass device	21 kg (46 lb)	23 kg (51 lb)
	Maximum rated load	53 kg (117 lb)	70 kg (154 lb)

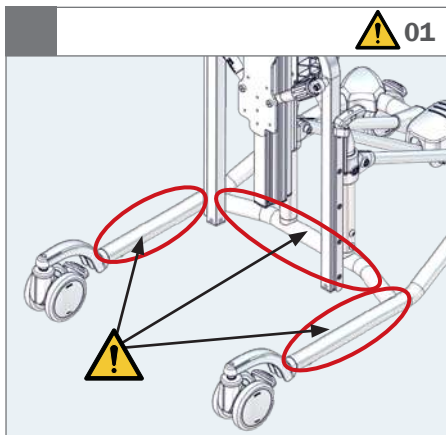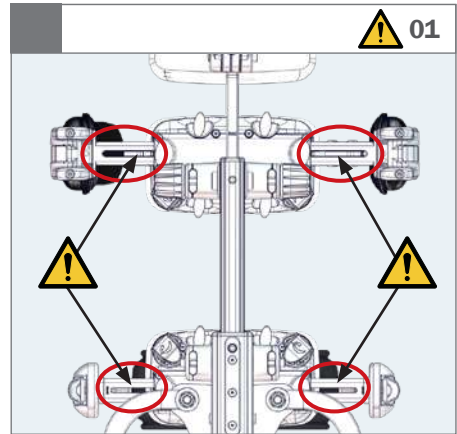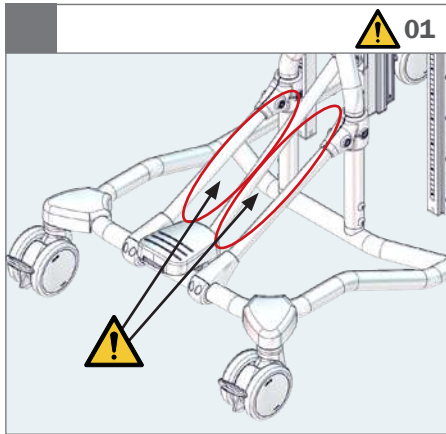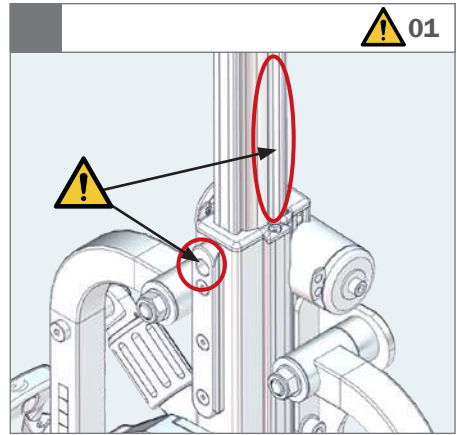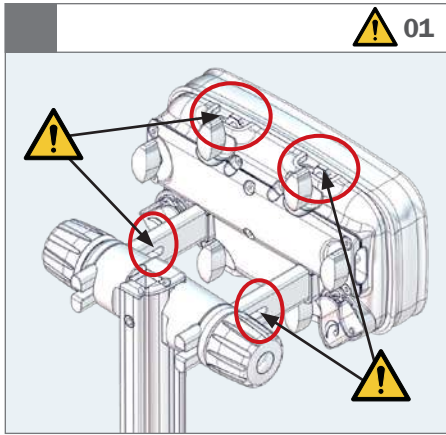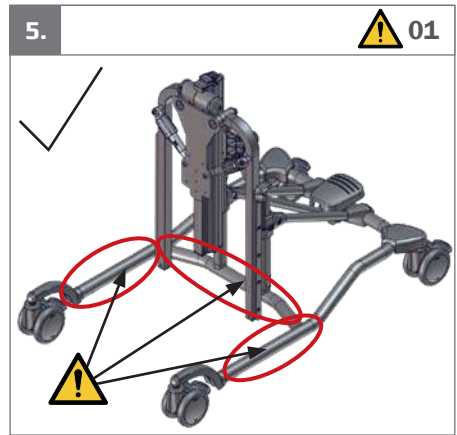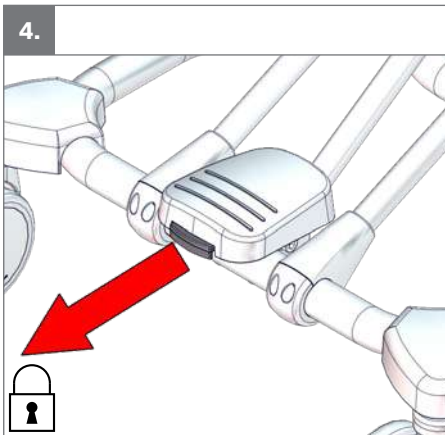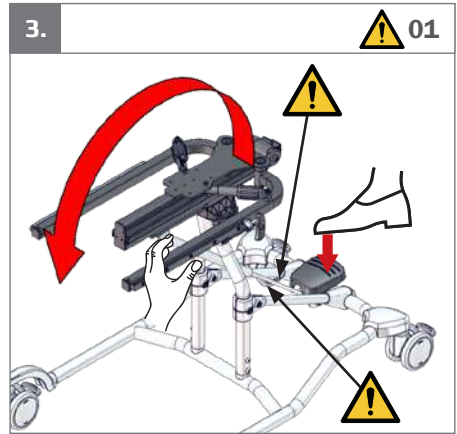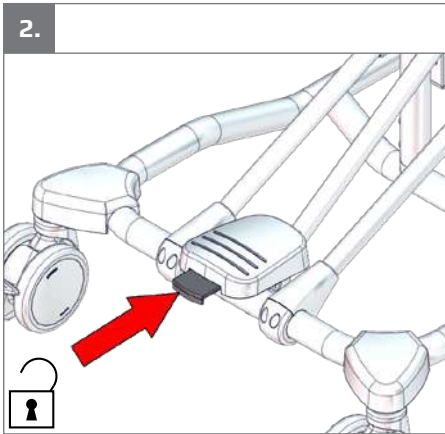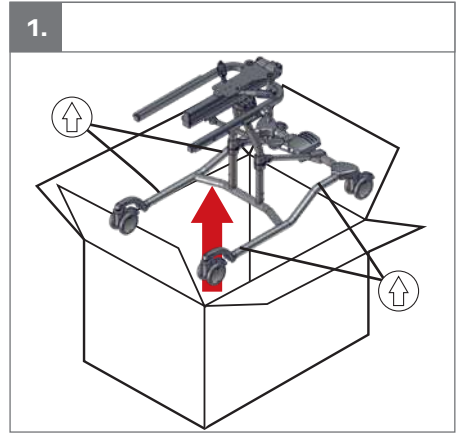
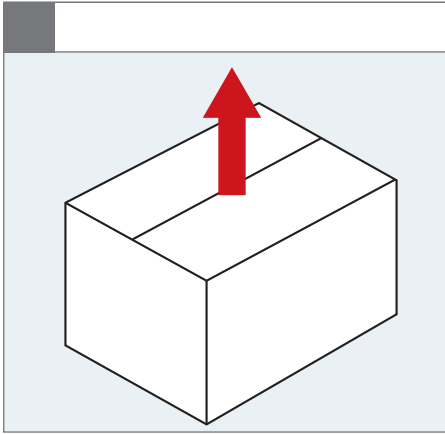
Flat packed

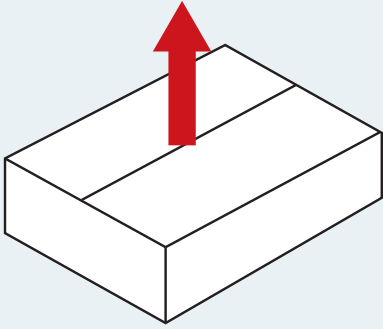
1.

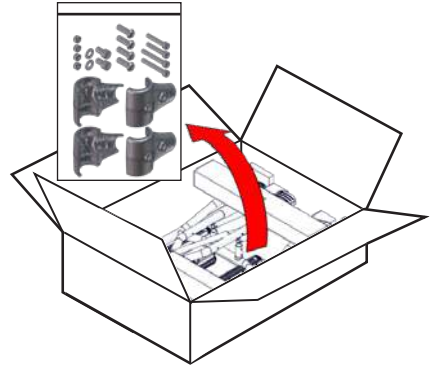
2.

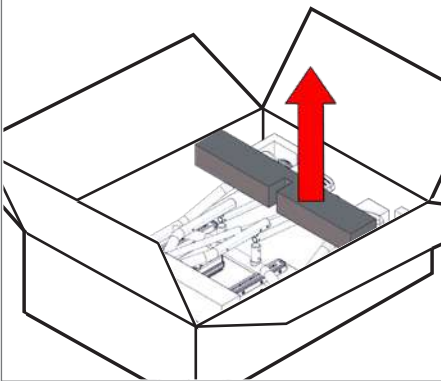
3.












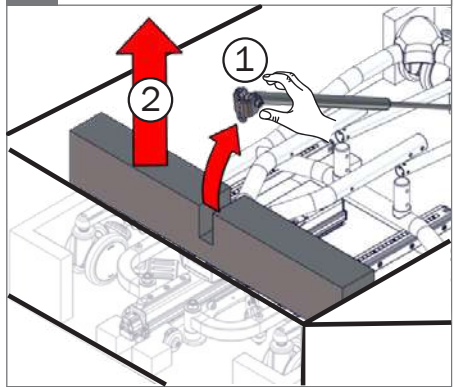
4.




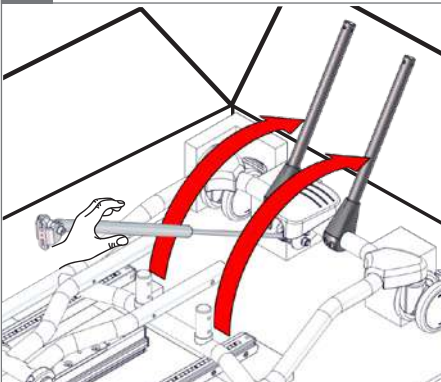
5.

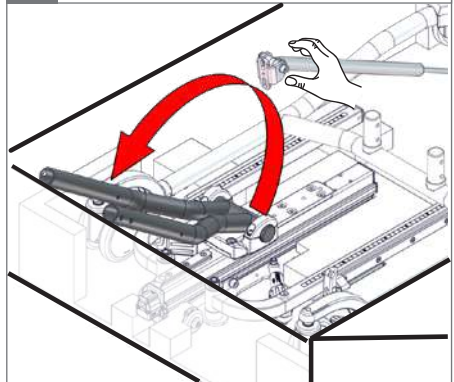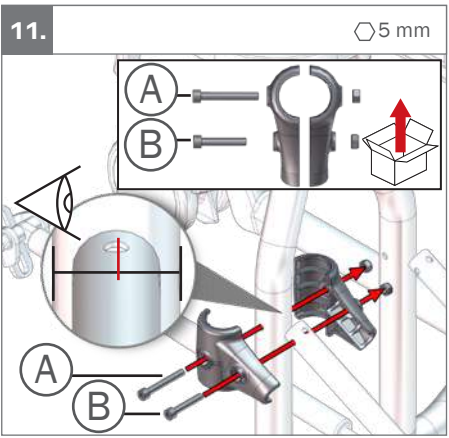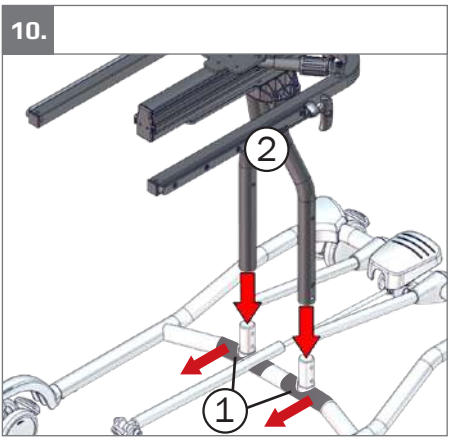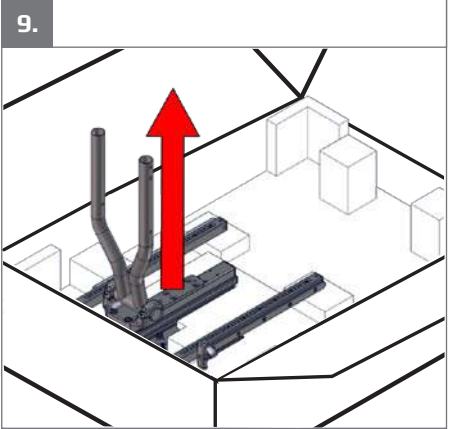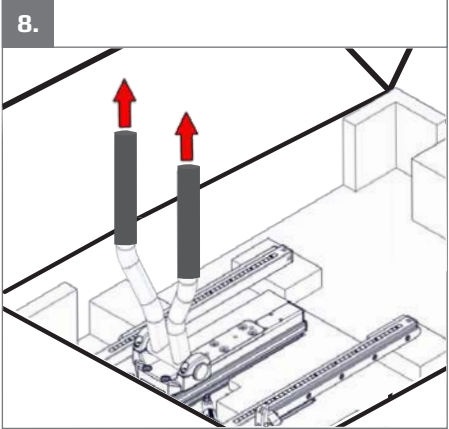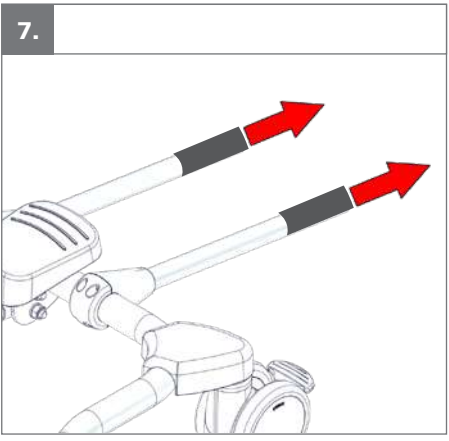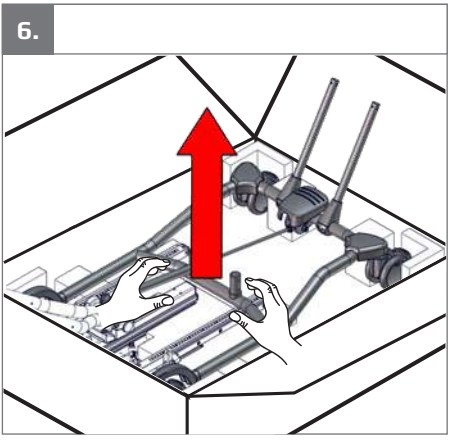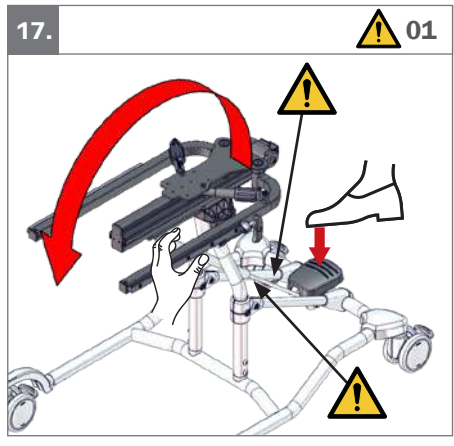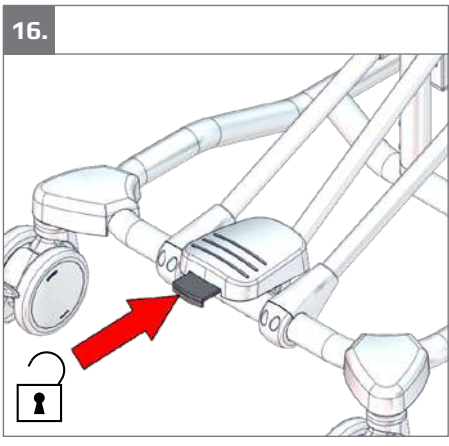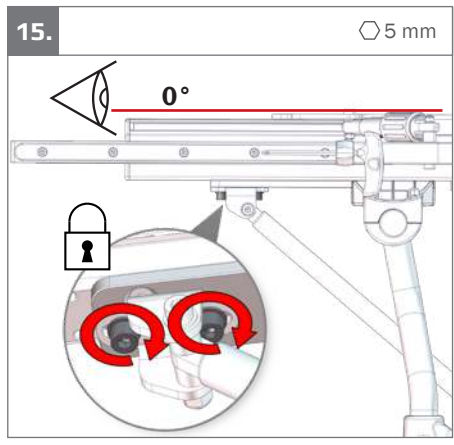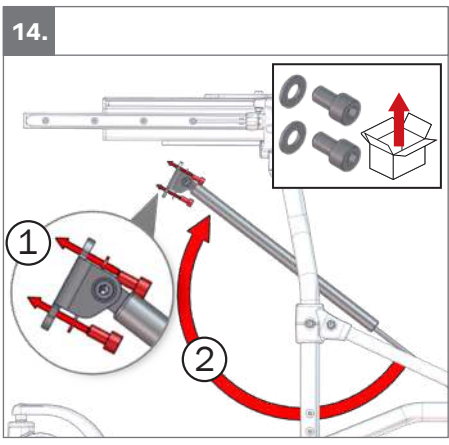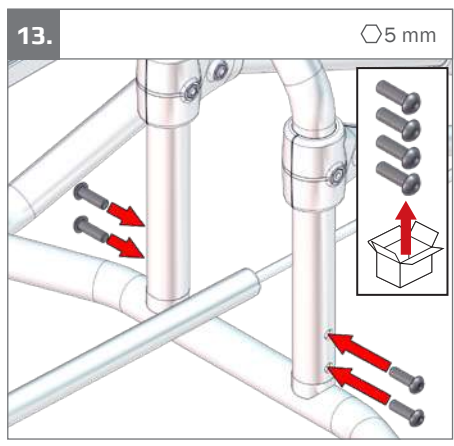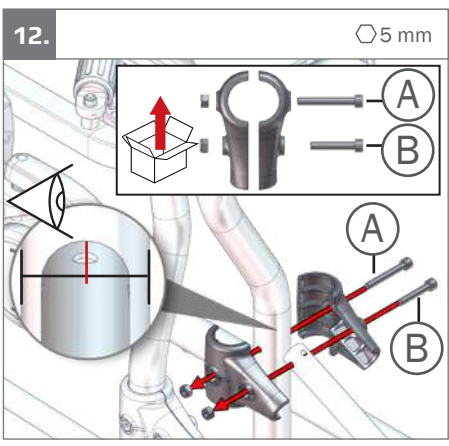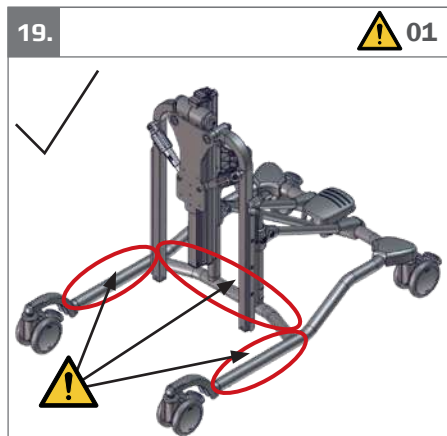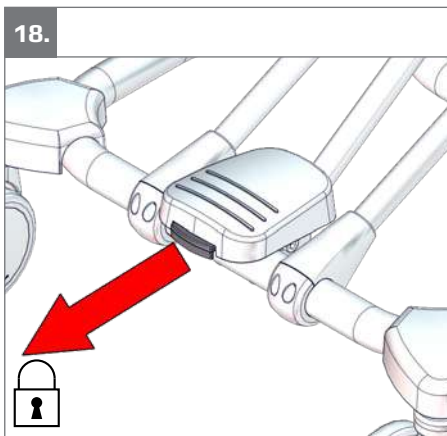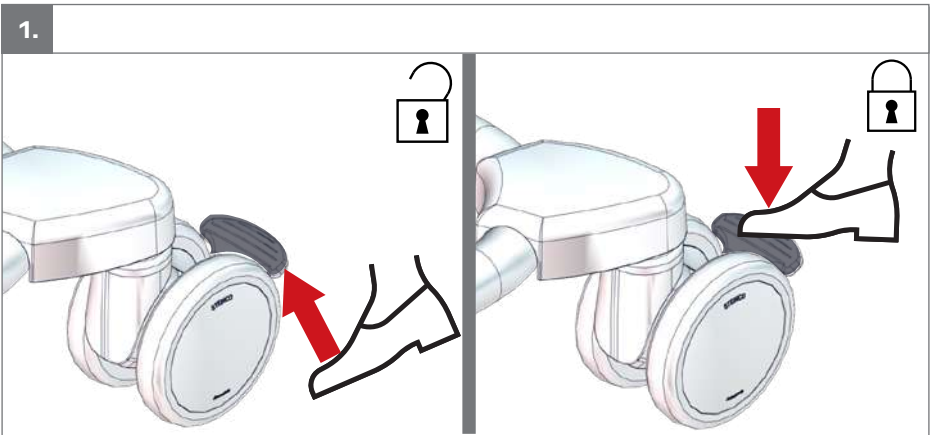
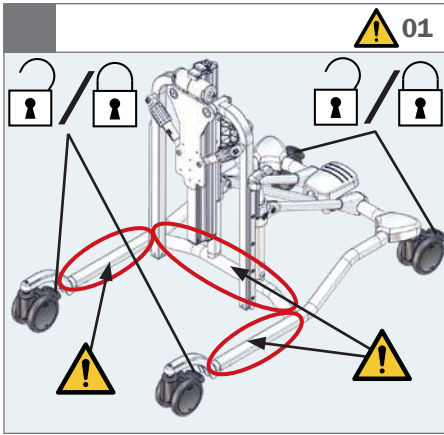
⚠ 01

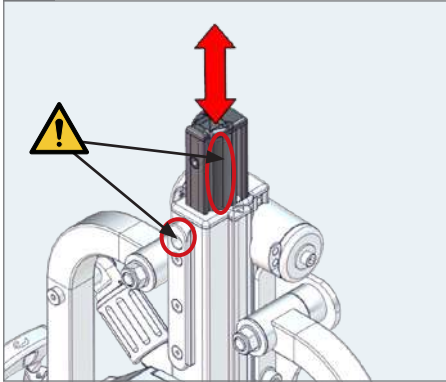1.



⊙ 5 mm












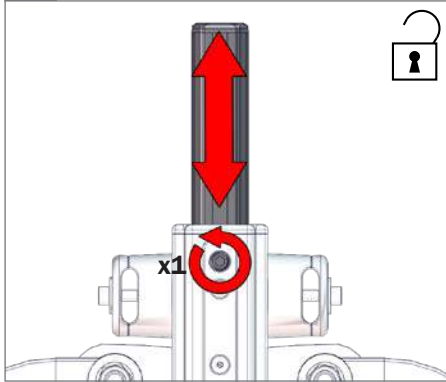
2.







⊙ 5 mm

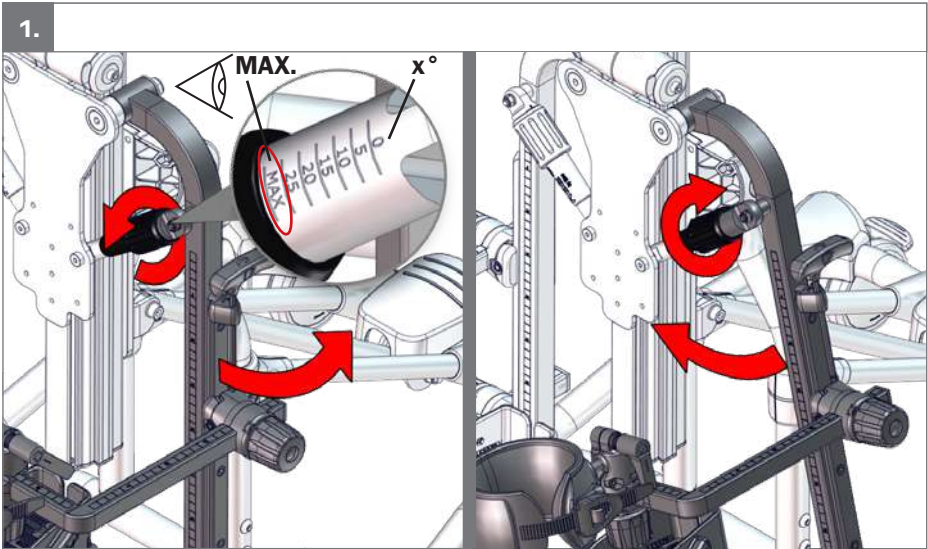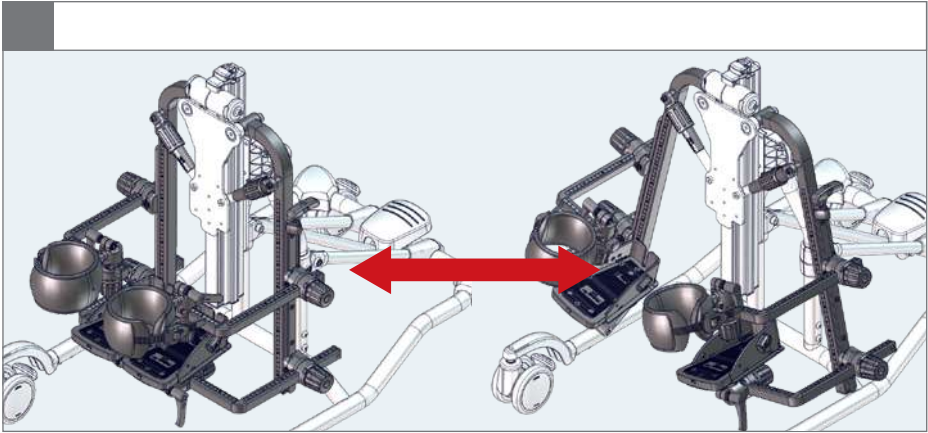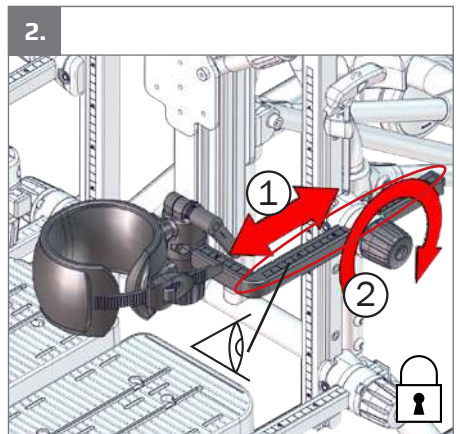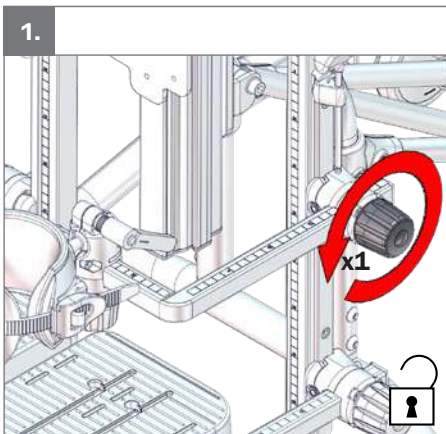
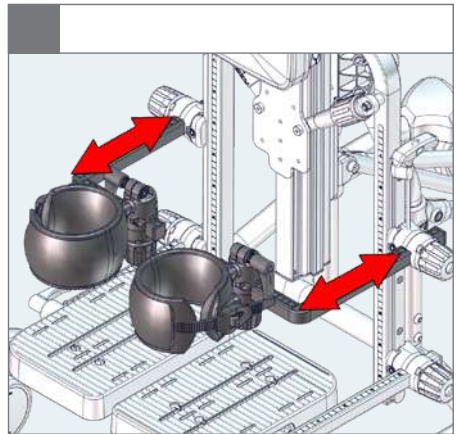
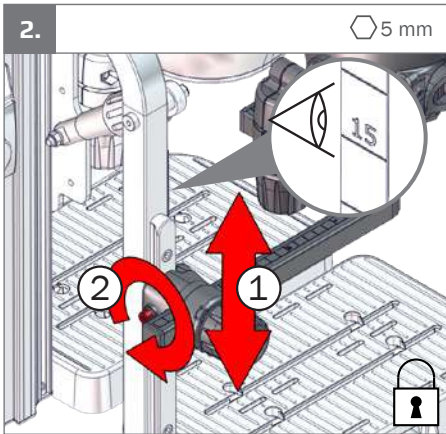
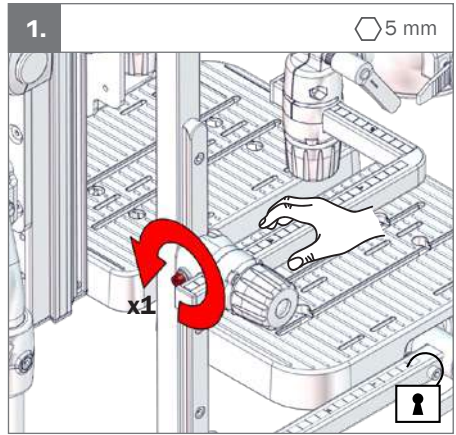
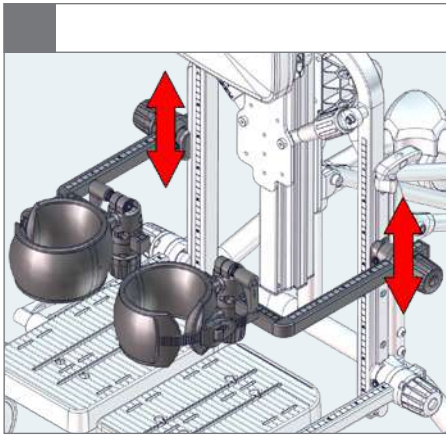
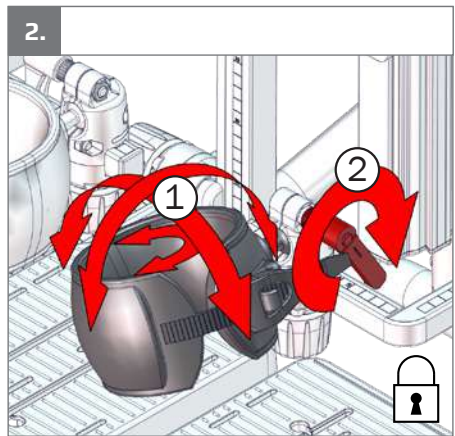
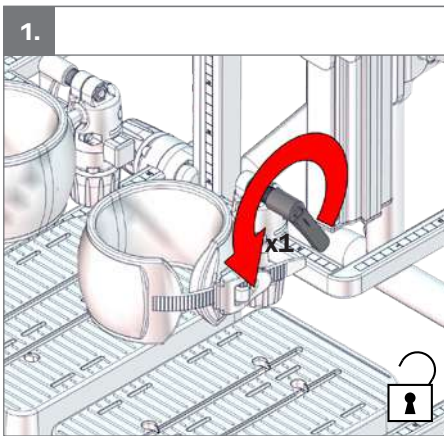
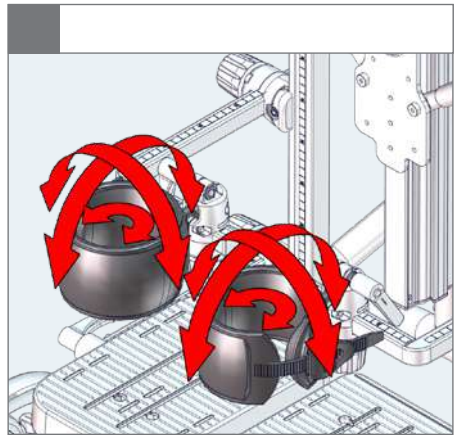
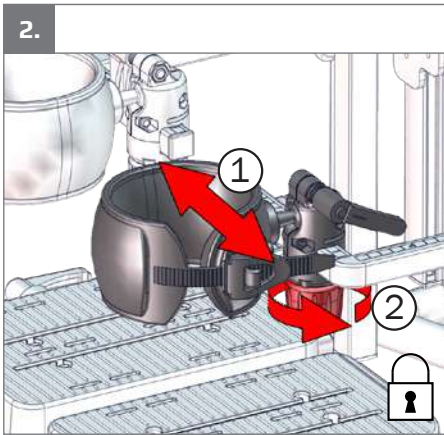
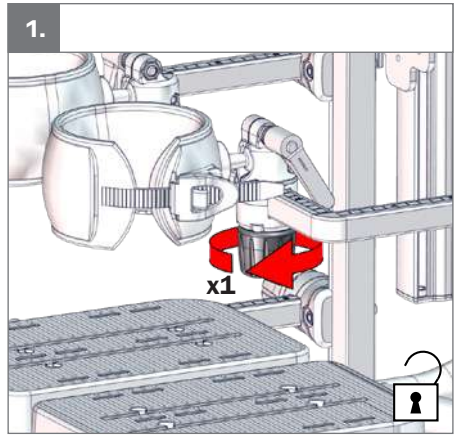
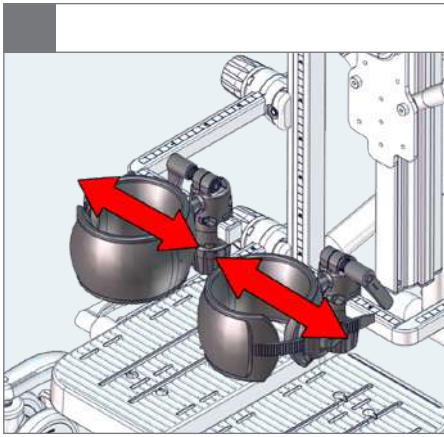
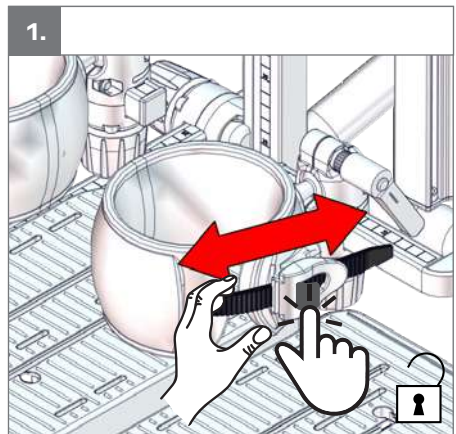
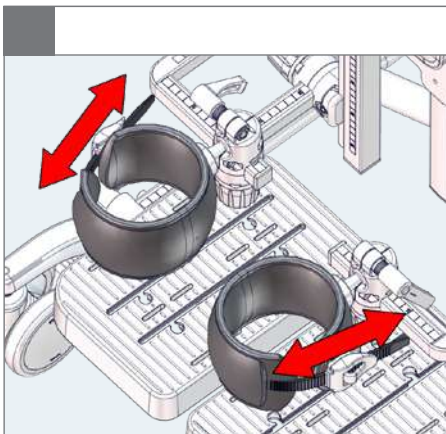
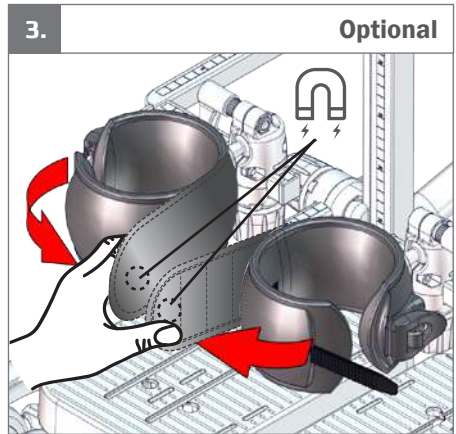
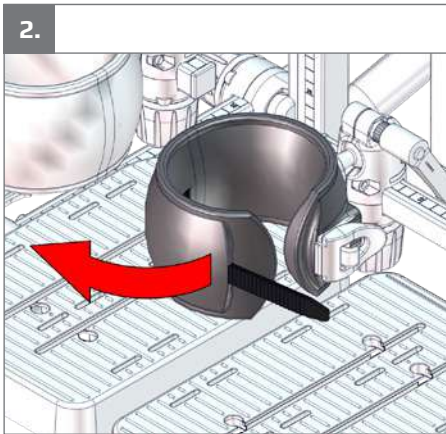
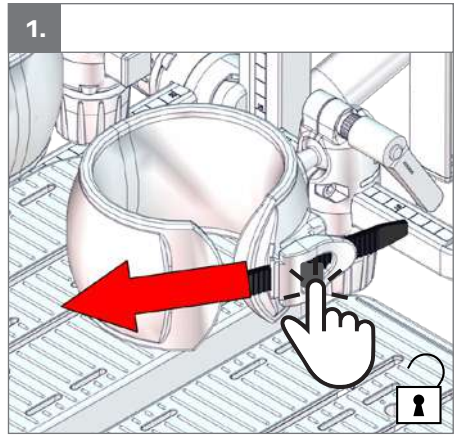
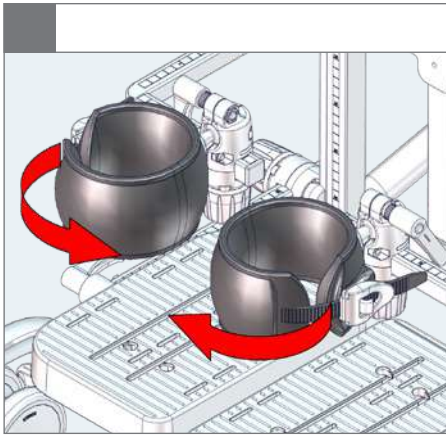
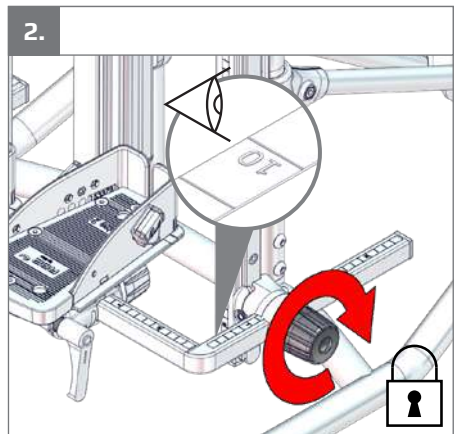
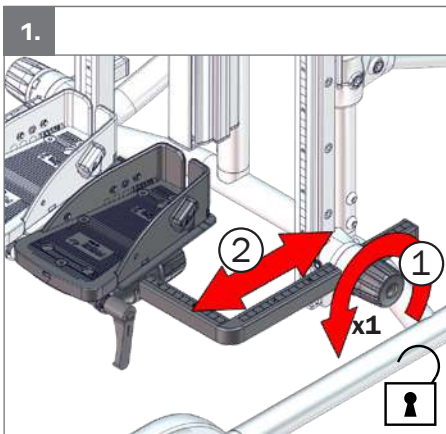
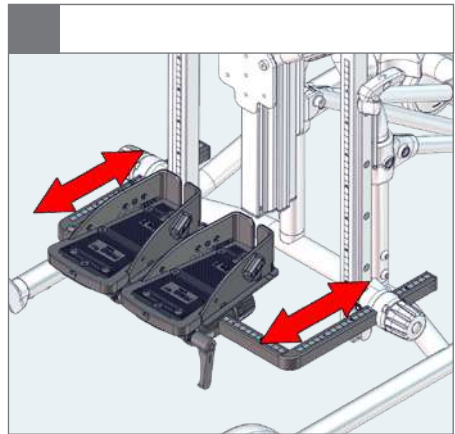
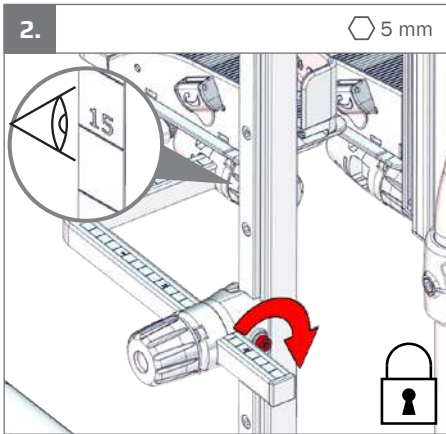
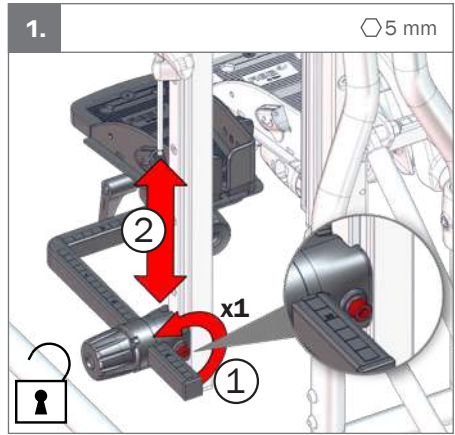
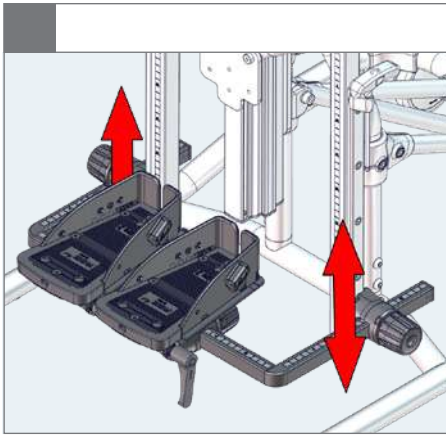
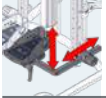
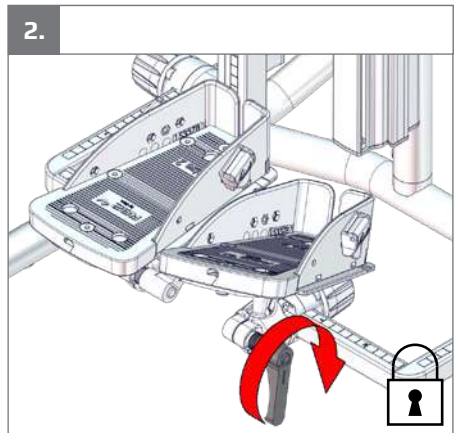
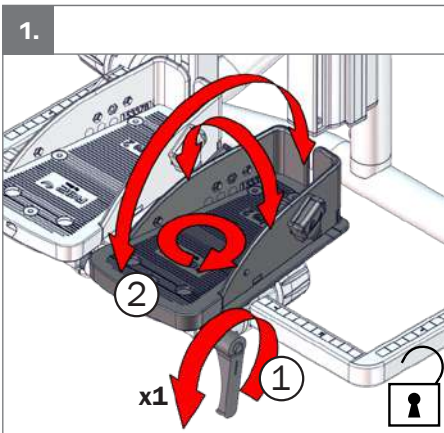
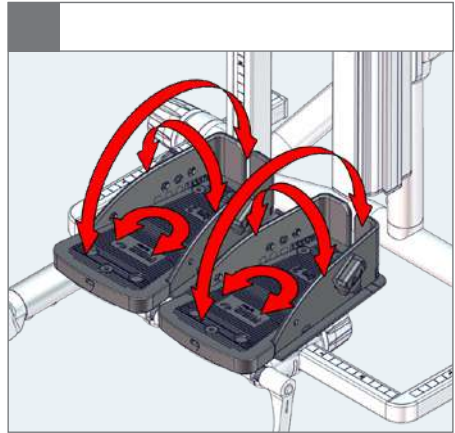
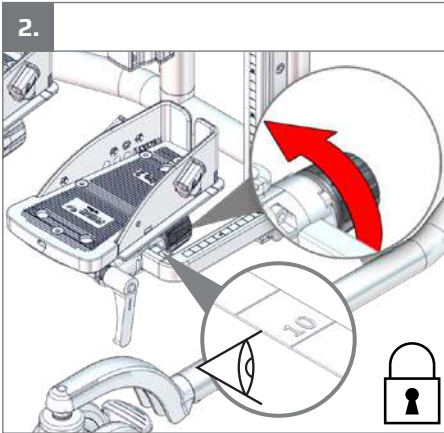
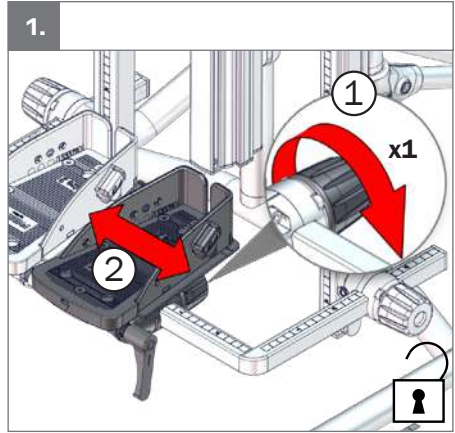
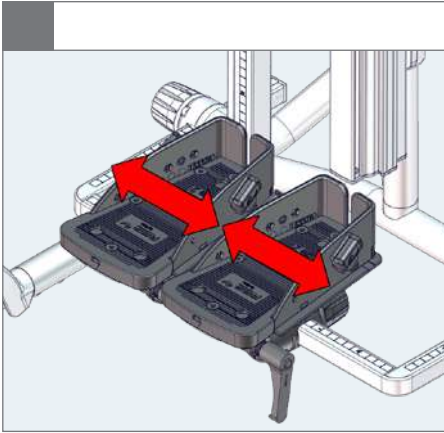
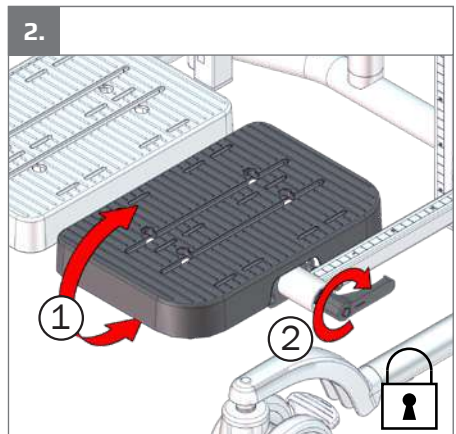
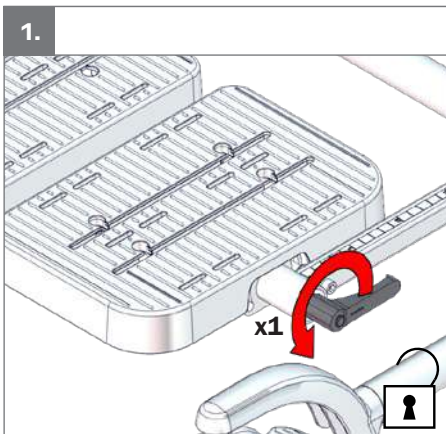
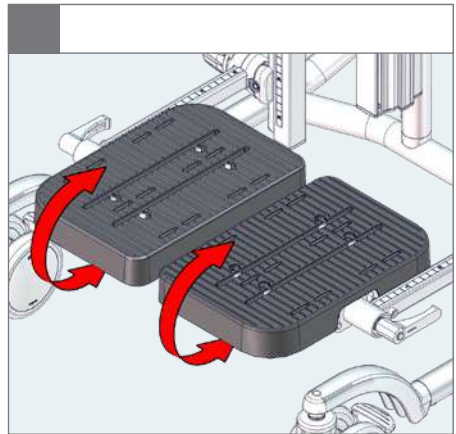
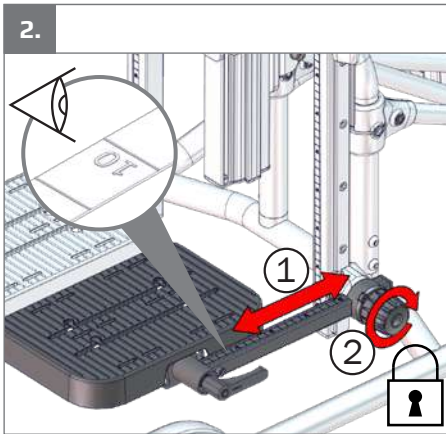
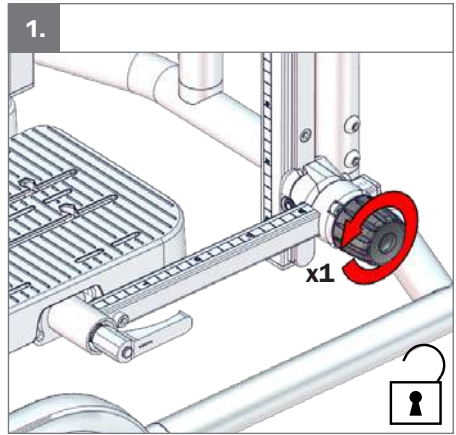
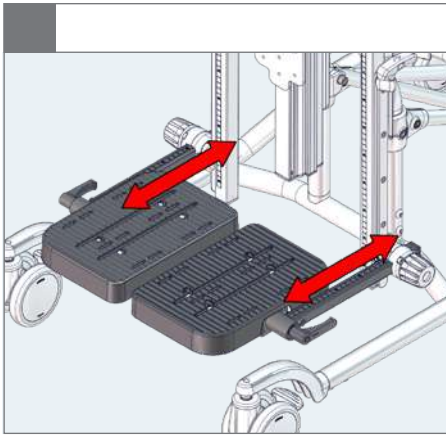
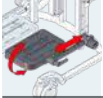
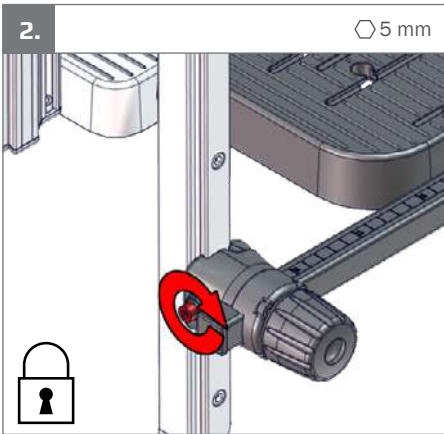
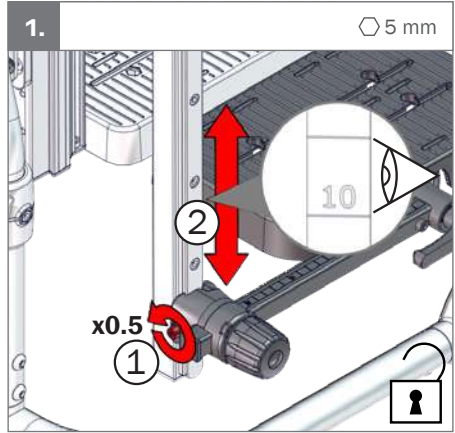
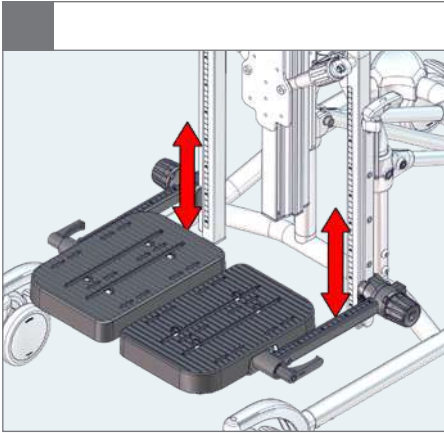
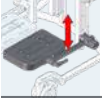
! 01

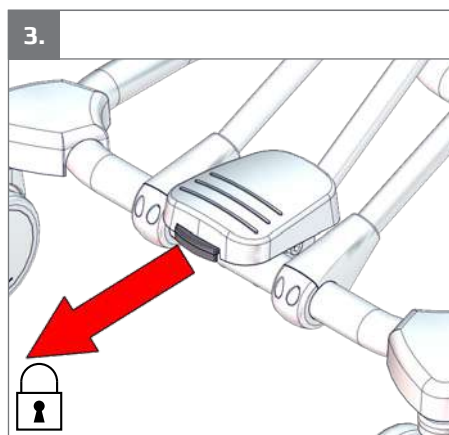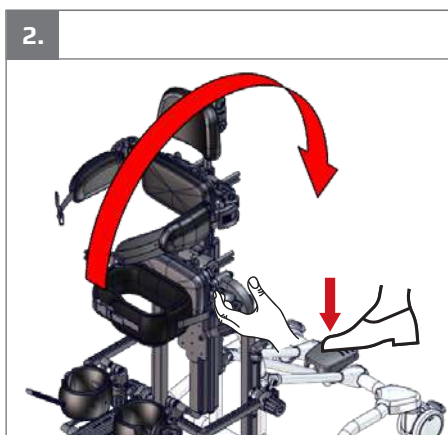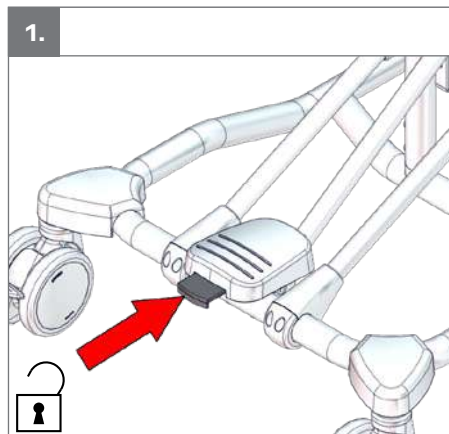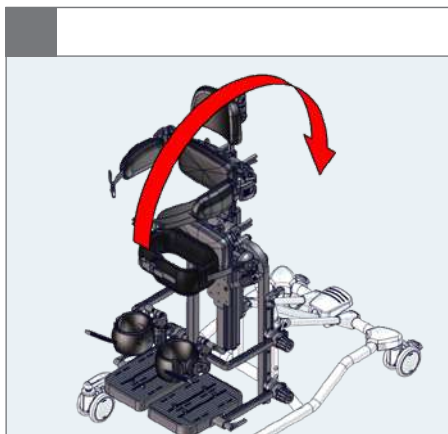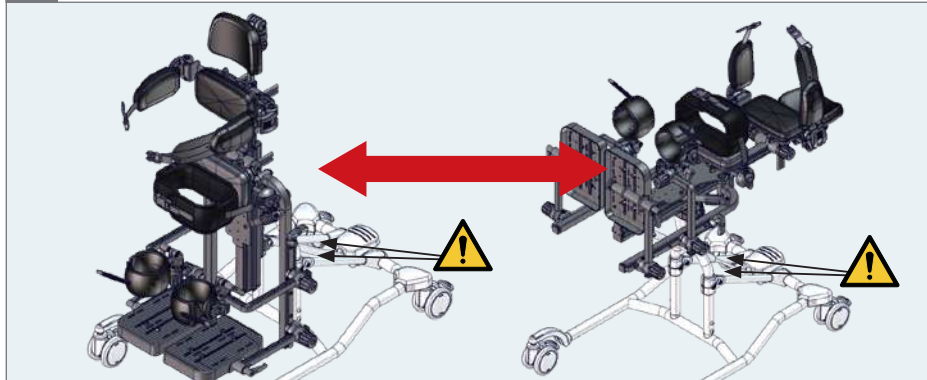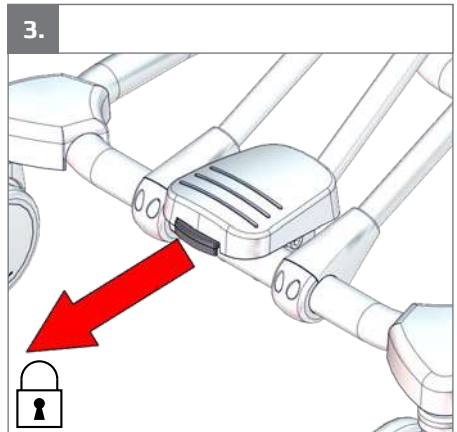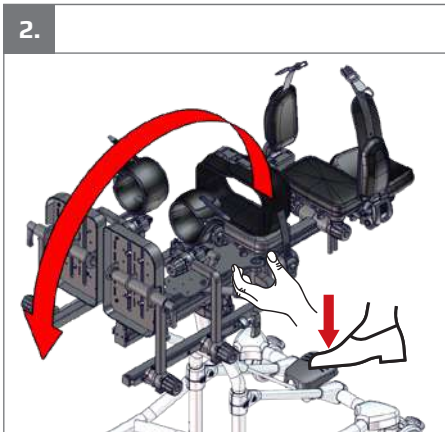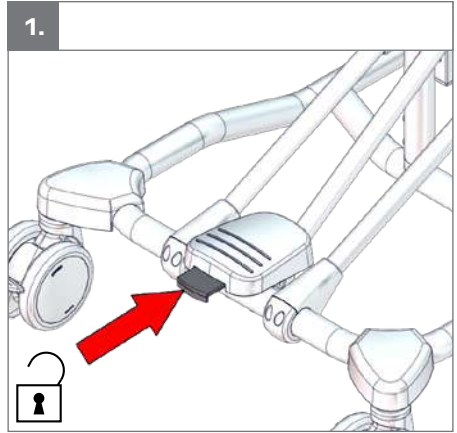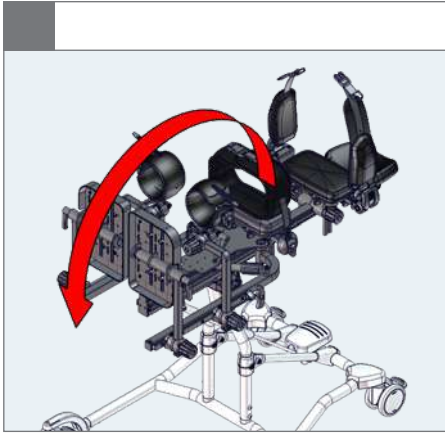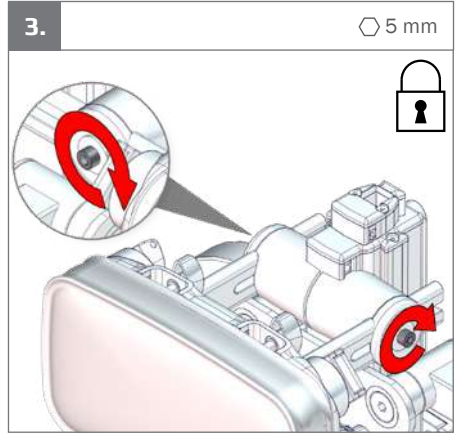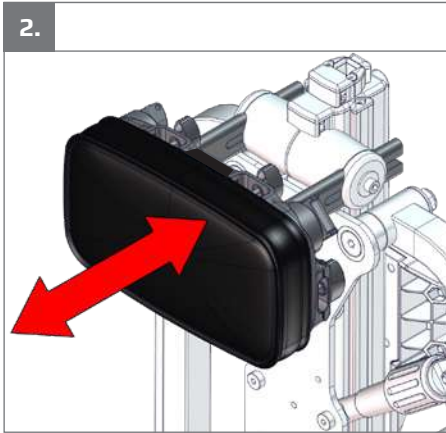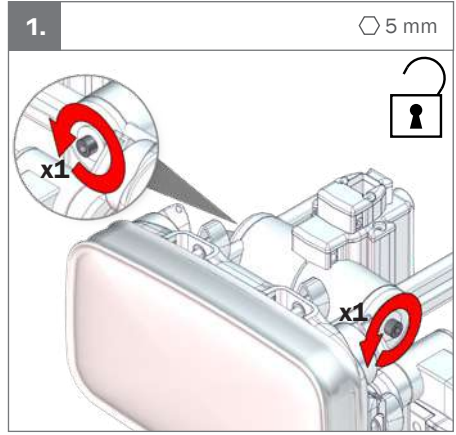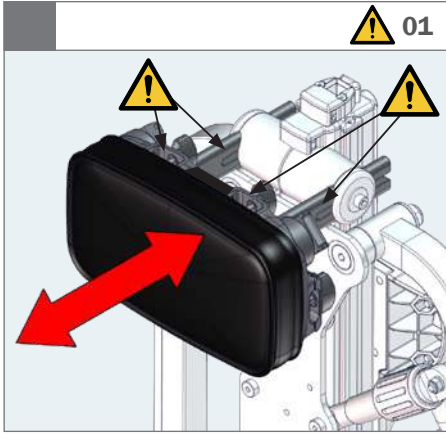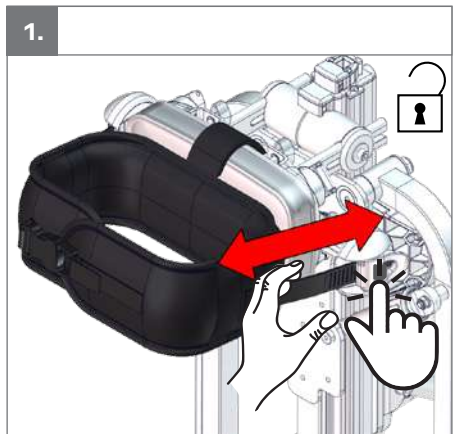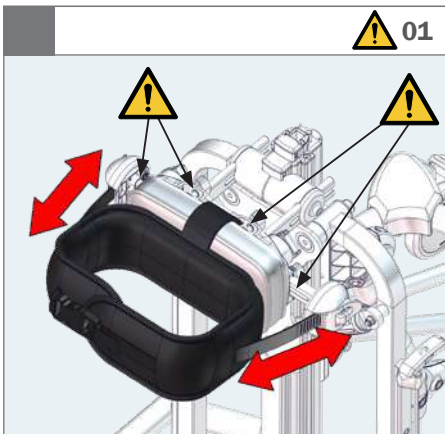
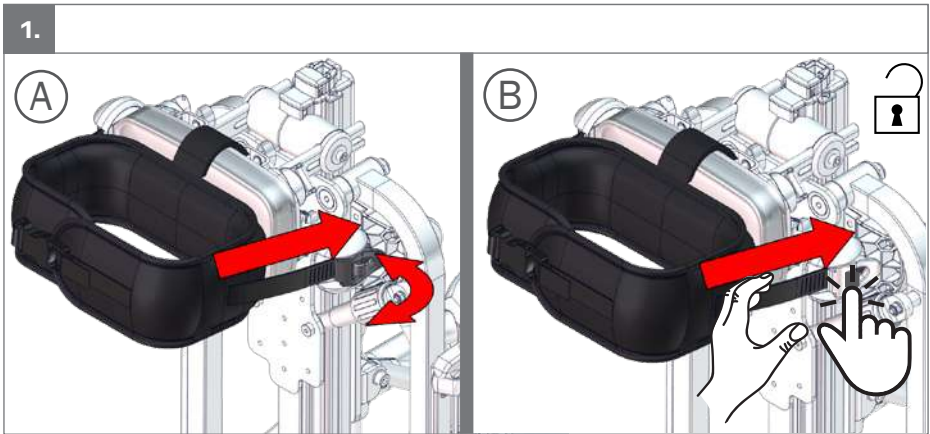
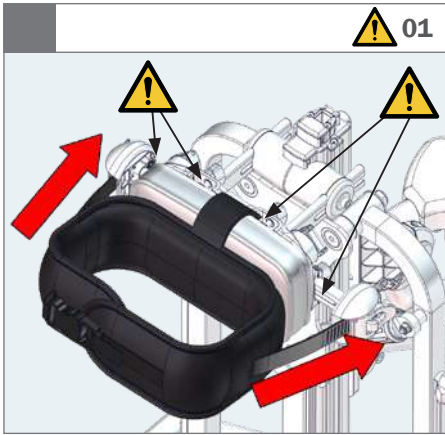
! 01

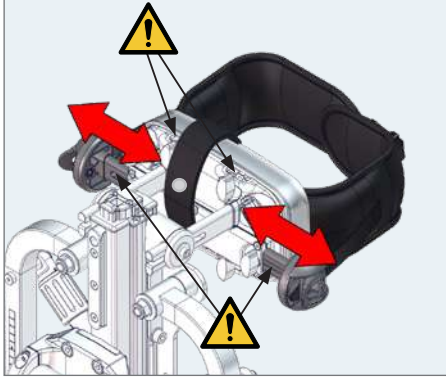
1.






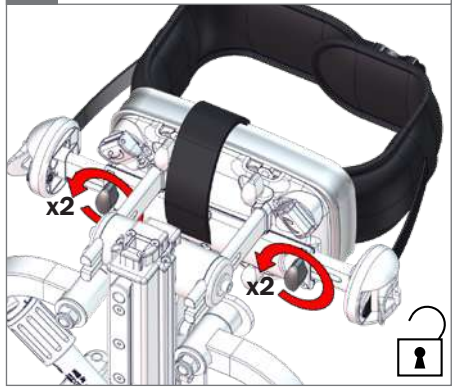
2.








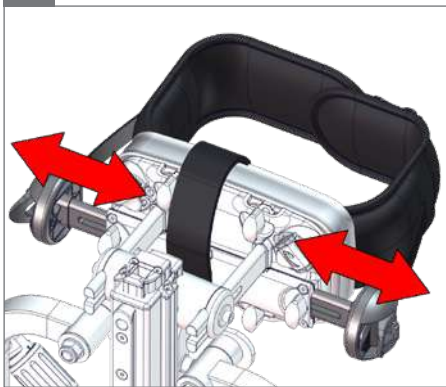
3.













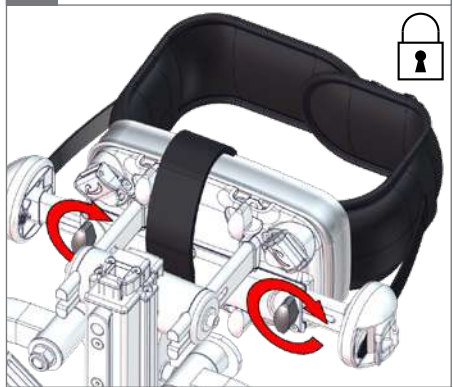
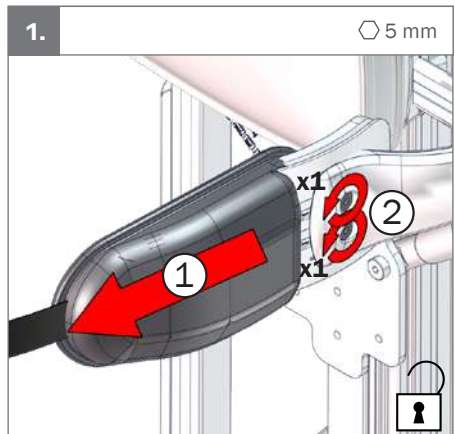
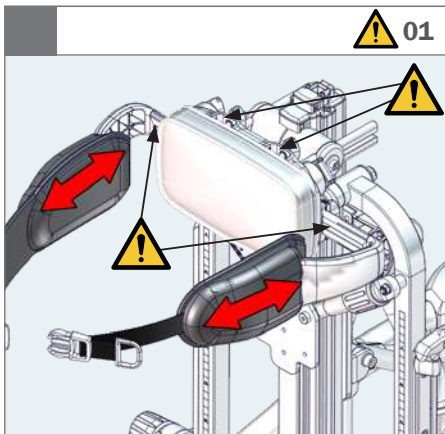
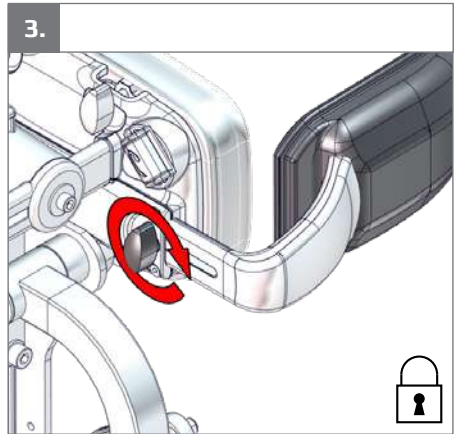
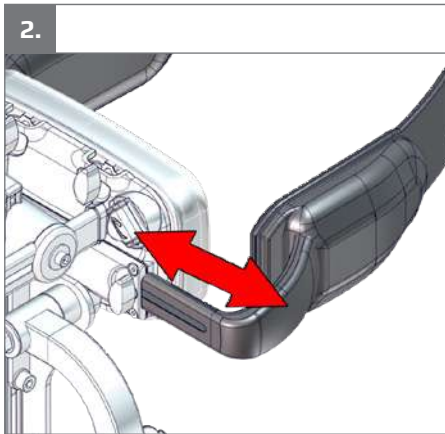
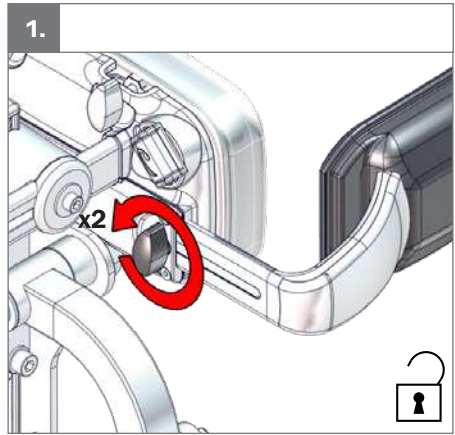
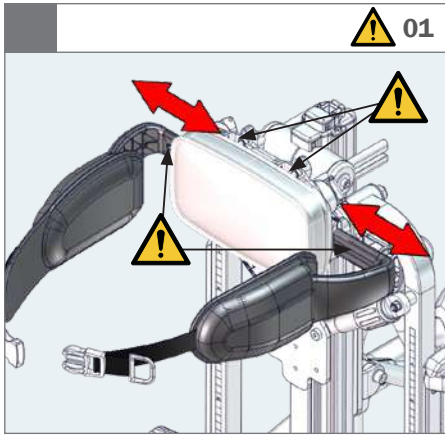
5 mm

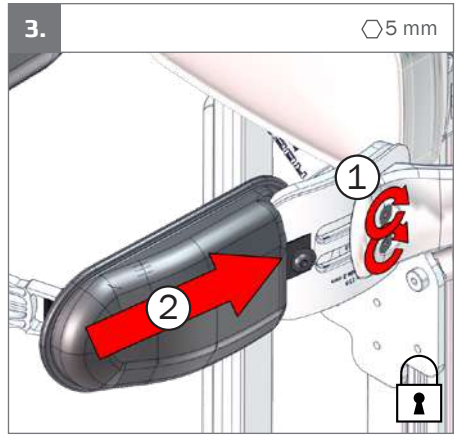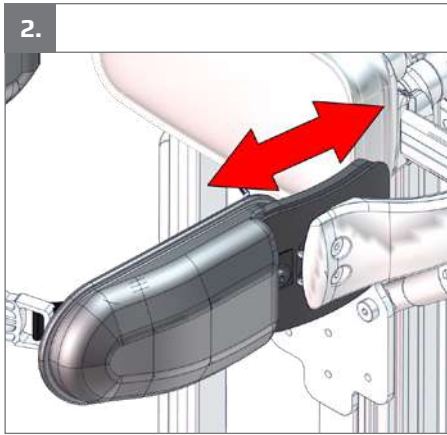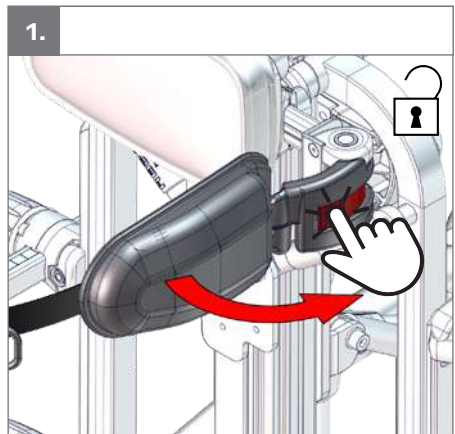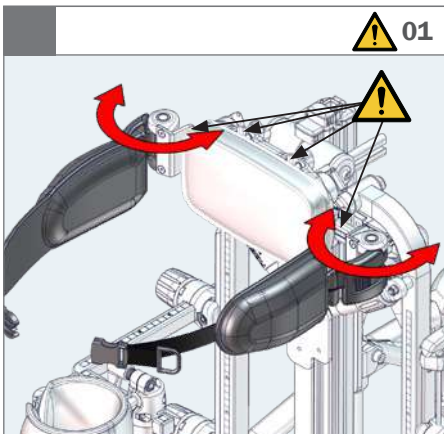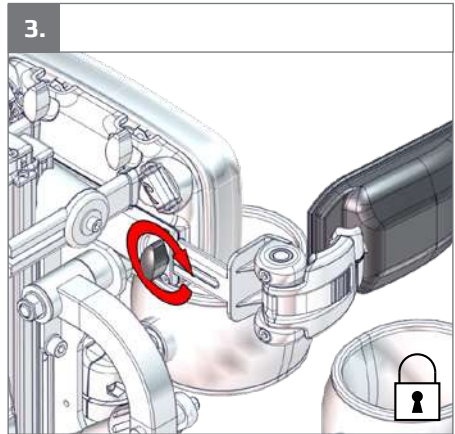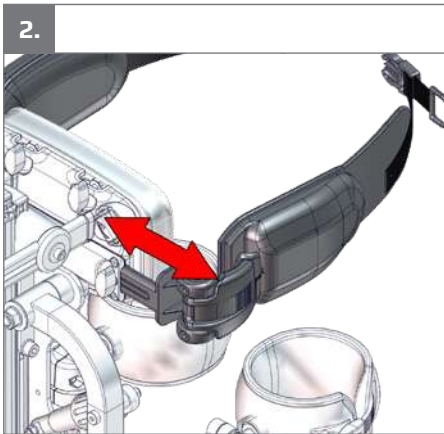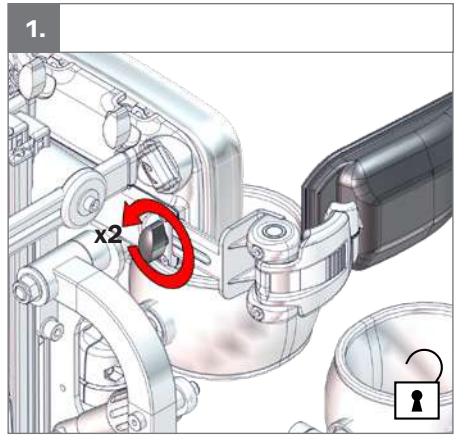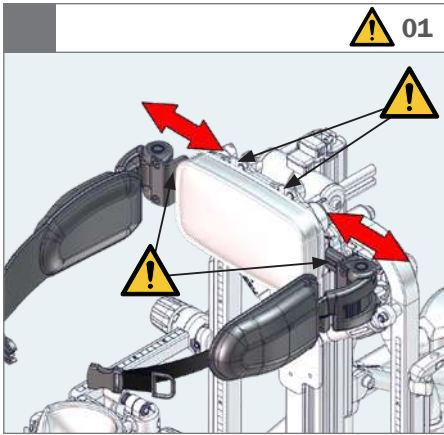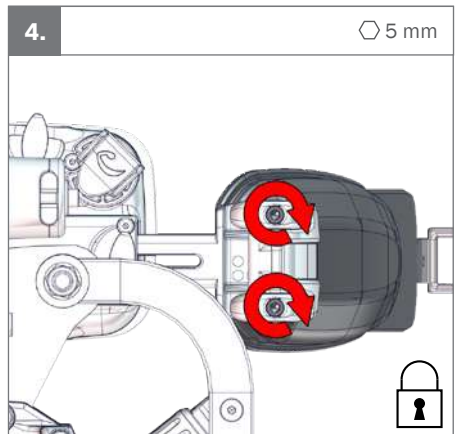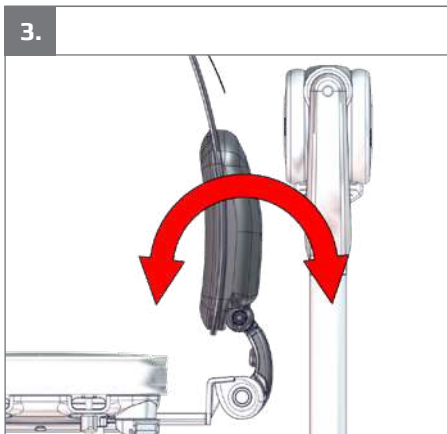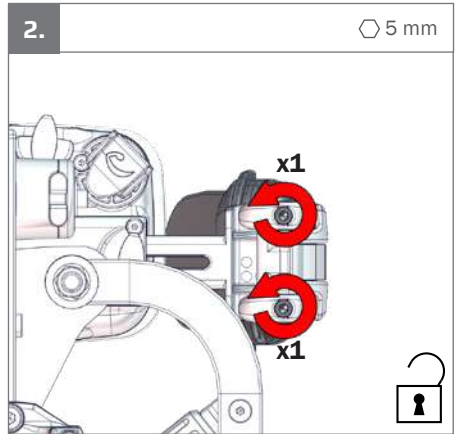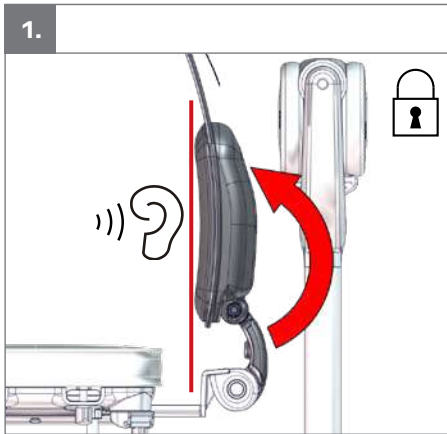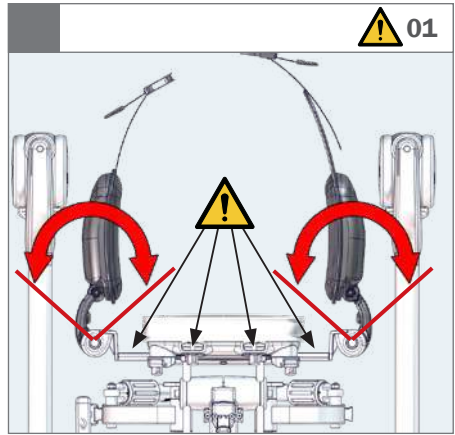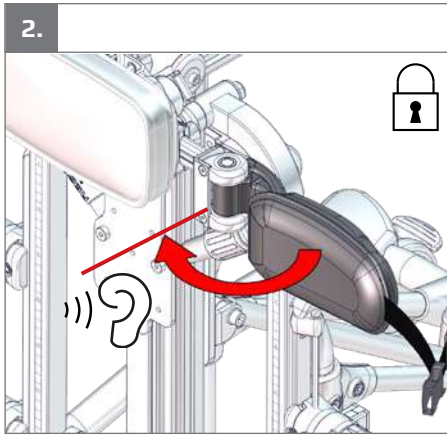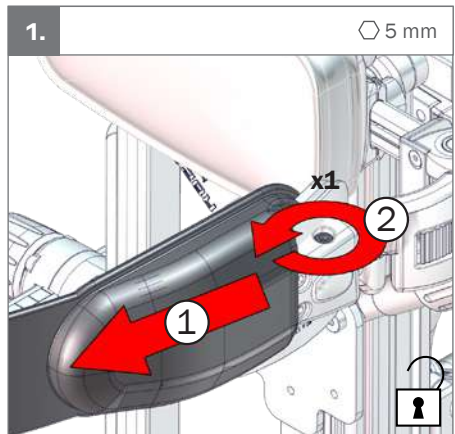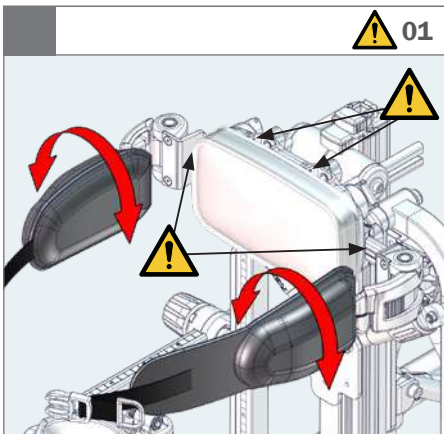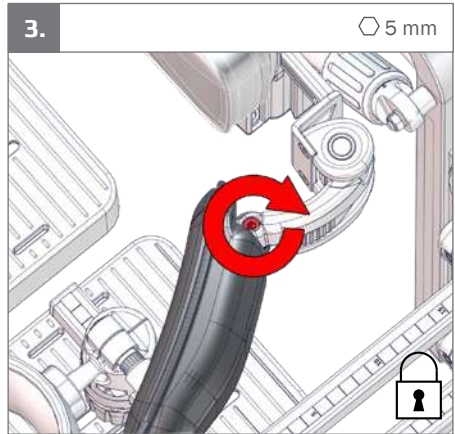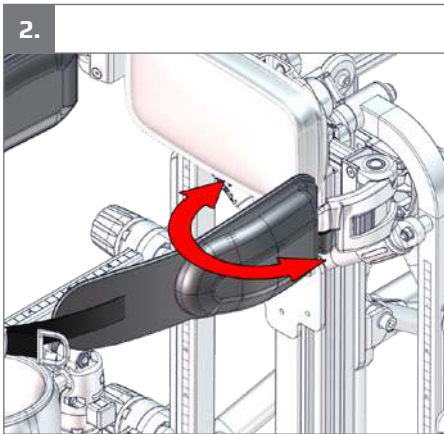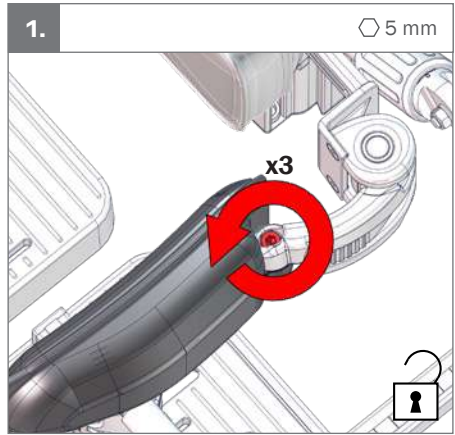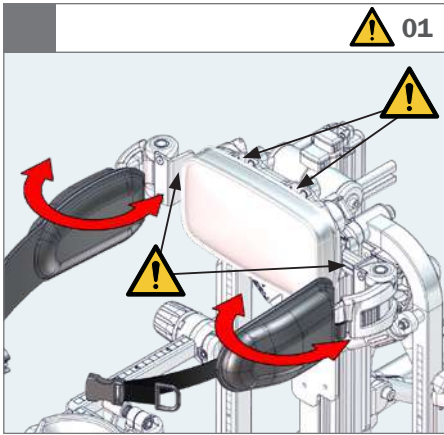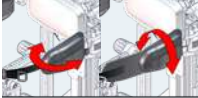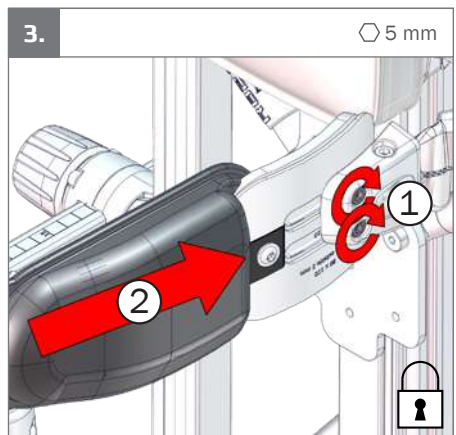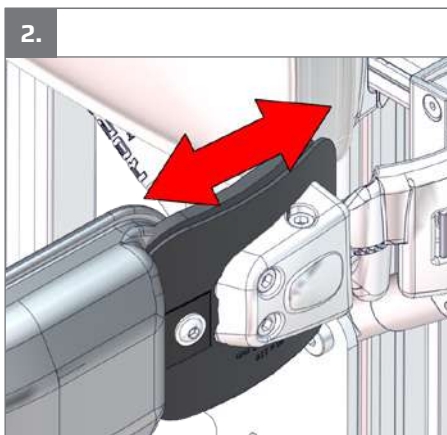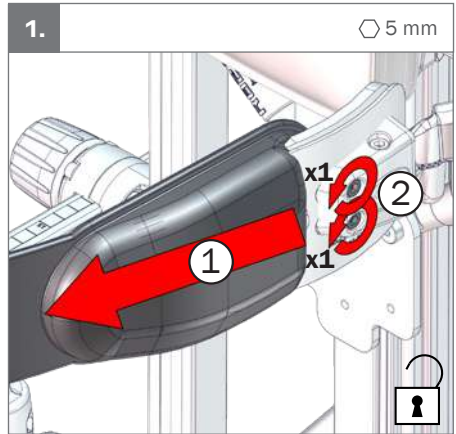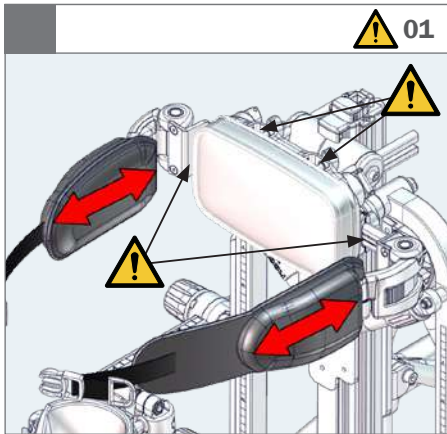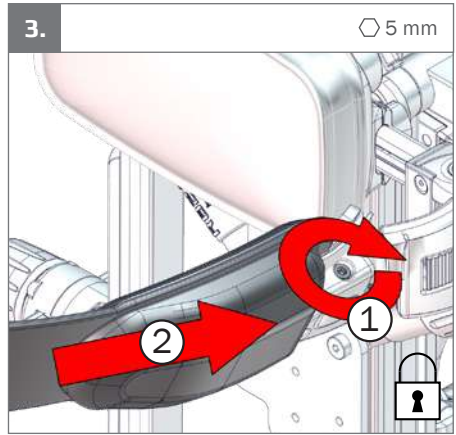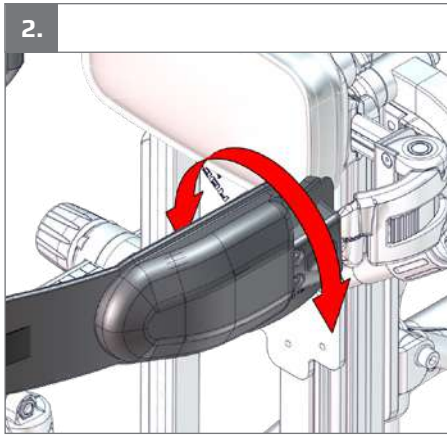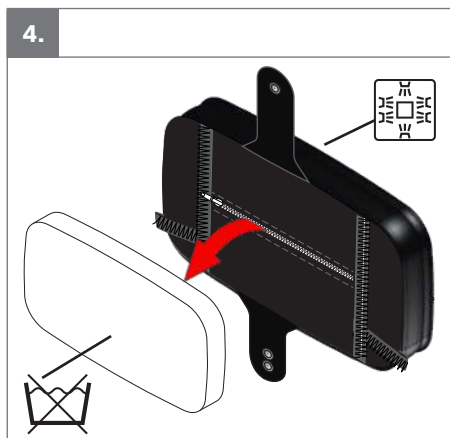
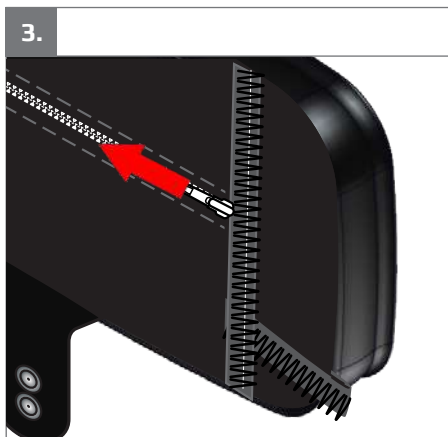
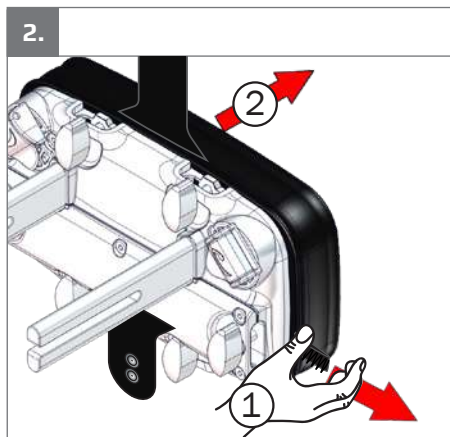
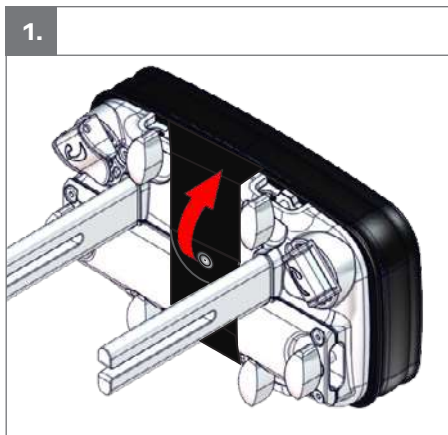
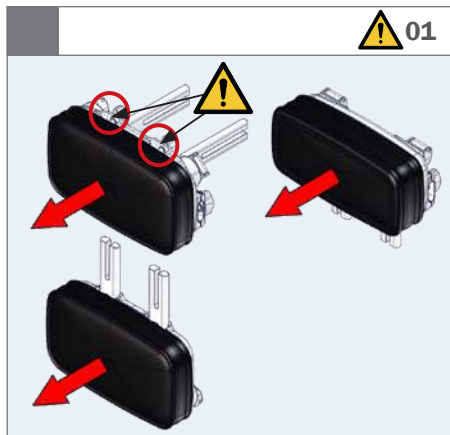
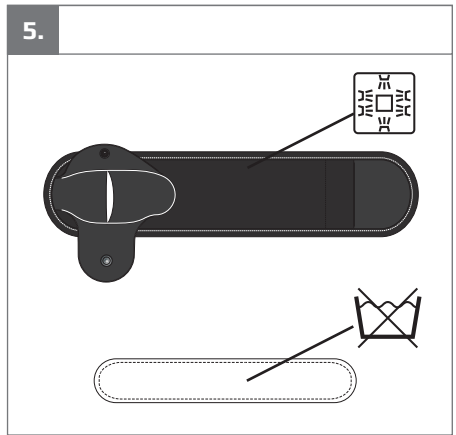
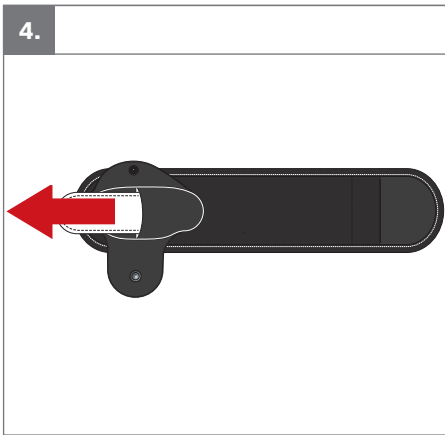
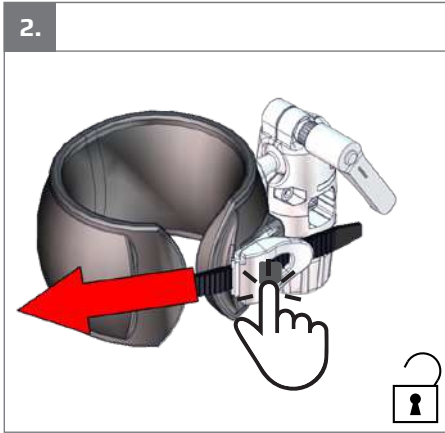
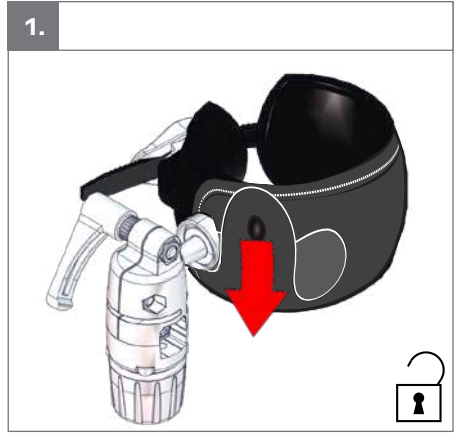
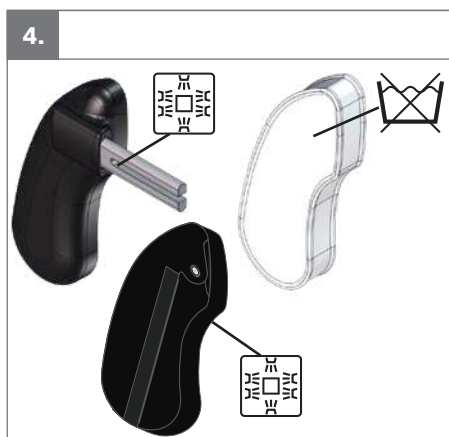
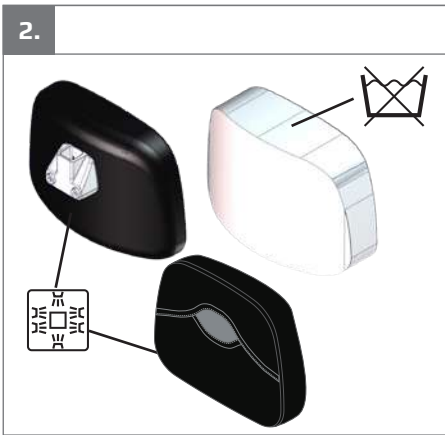
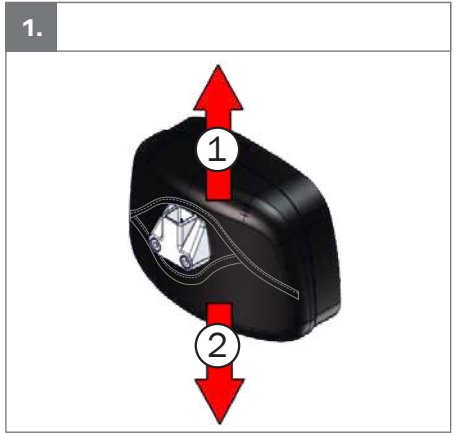
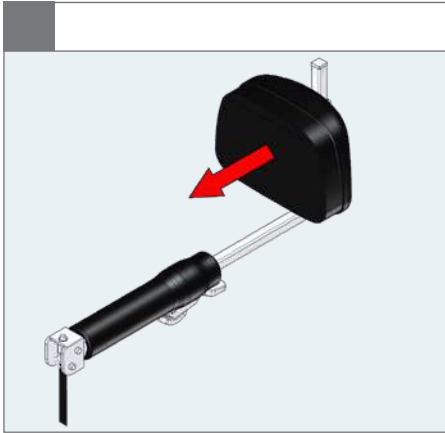
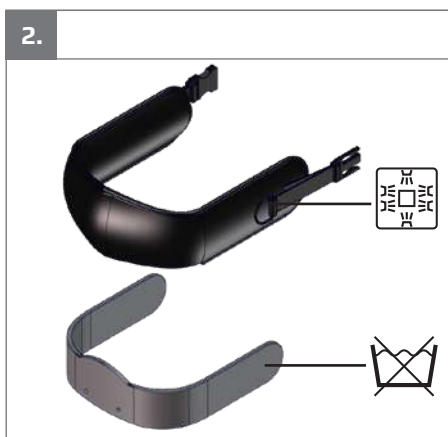
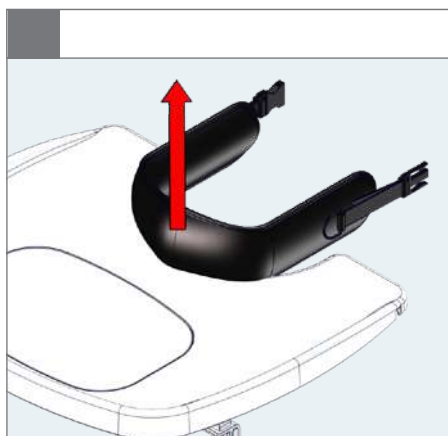
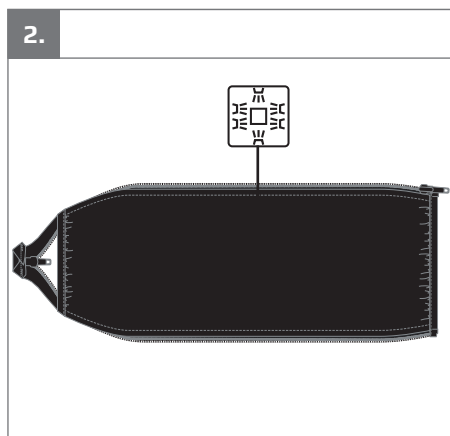
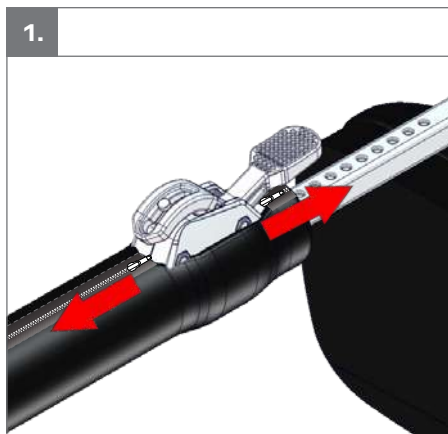
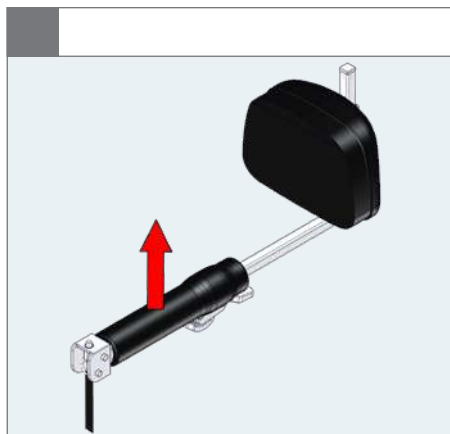
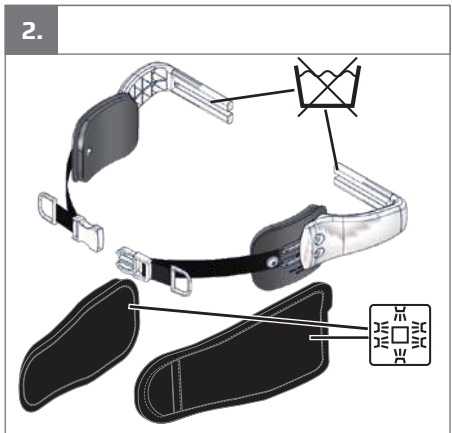
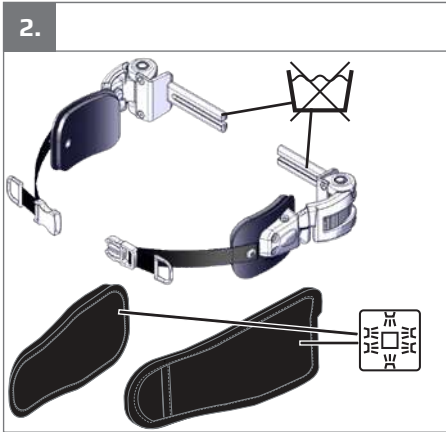
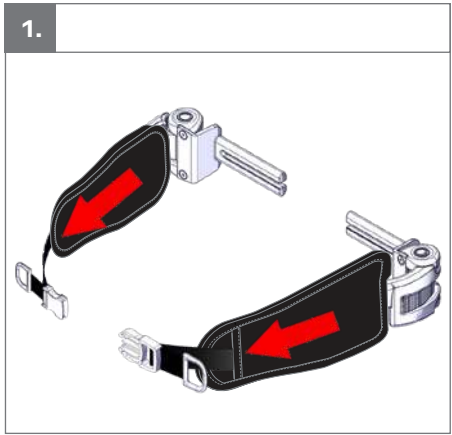
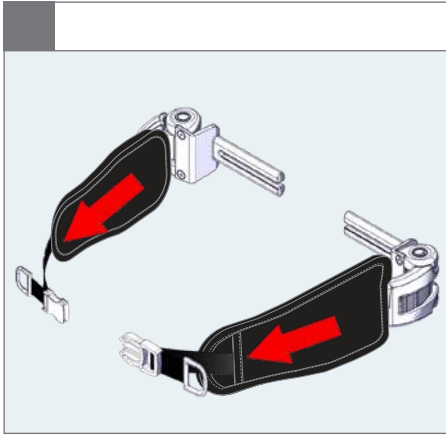
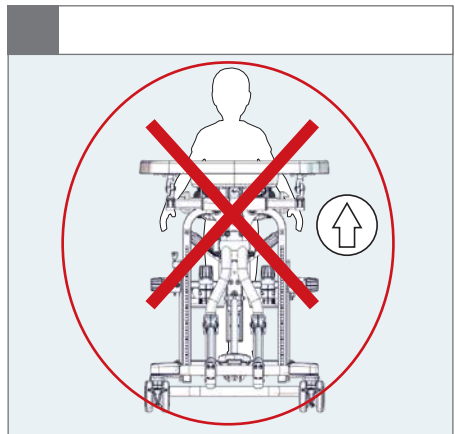
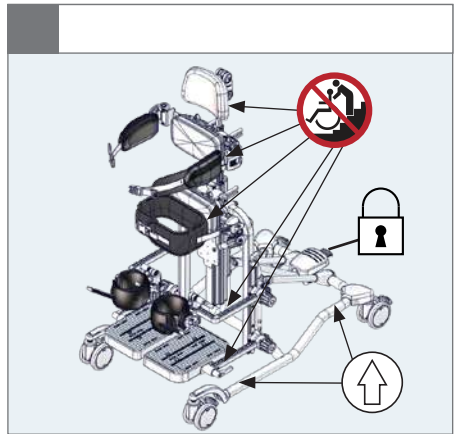
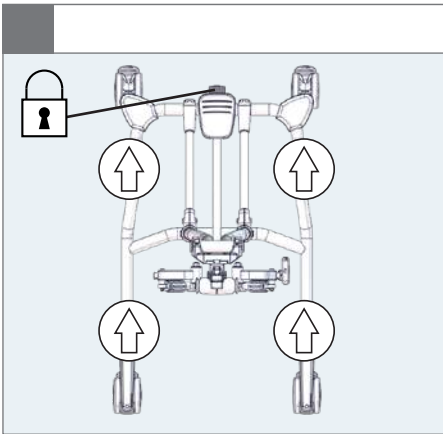
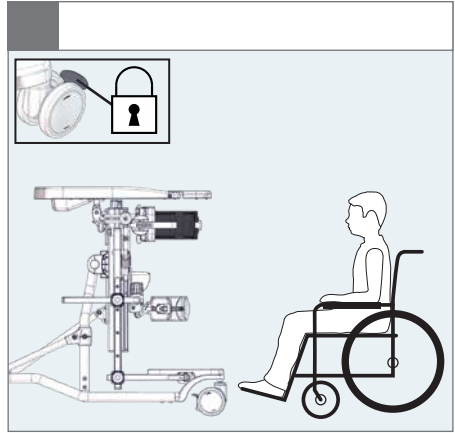
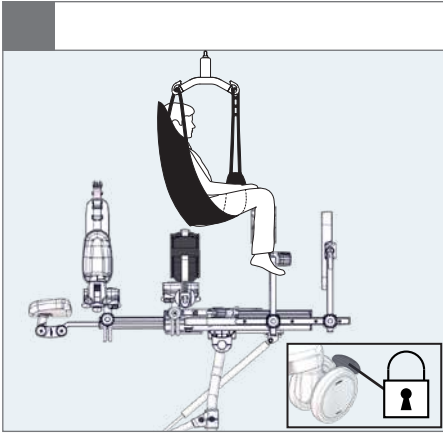
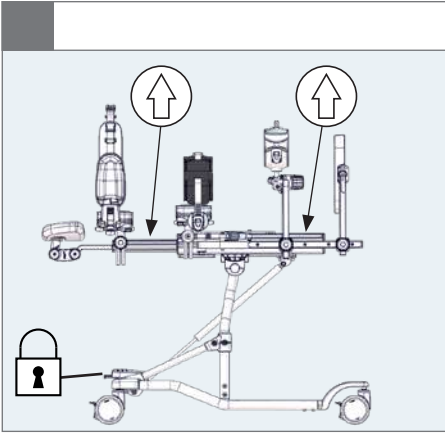
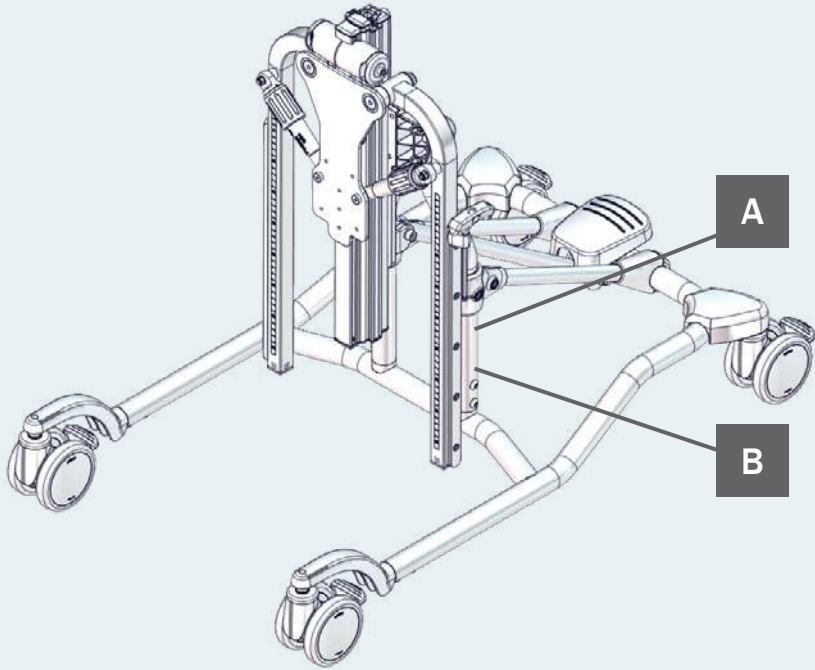
**A**

**R82**  
 by Etac

R82 A/S, Parallelvej 3  
 DK-8751 Gedved  
 etac.com  
 Made in xxxxxxxx

MD

CE UK CA

YYYY-MM-DD

**B**

Posting group  
 Product name, size x


DMR ver.: xxxx

SW xxxxxxxxxxxx


REF xxxxxxxxxxxxxxxxxxxx

▲ Max. user mass: xxx kg  
 ▲ Max. rated load: xxx kg  
 ▲ Mass device: xxx kg


1001 (XX)XXXXXXXXXXXX(XX)XXXXXX(XX)XXXXXX

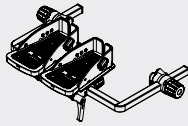
9996098259

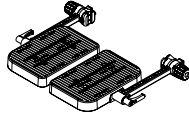
9996098258

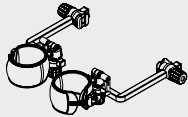
9996098260

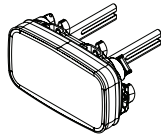
9996098257

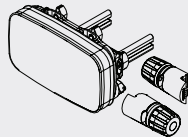9996098254

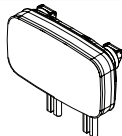
9996098255

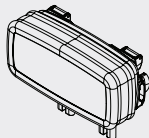
9996098256


9996098262

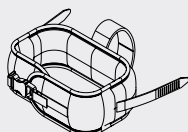
9996098261

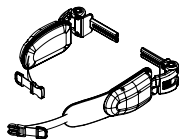
9996098263

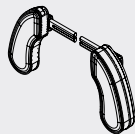
9996098270

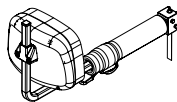
9996098264

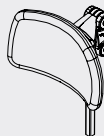
9996097080

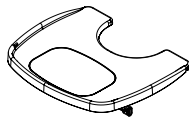
9996098265

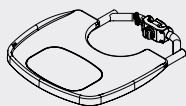
9996098269

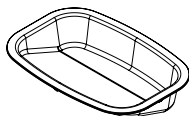
9996098268

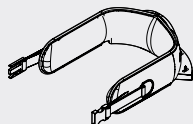
9996098266

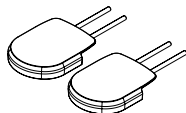
9996098267

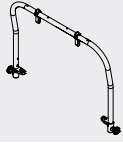
9996097706

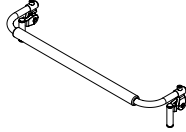
9996097717

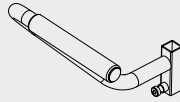
9996098271

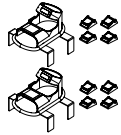
9996097241

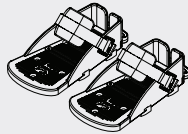9996097323

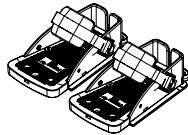9996098147

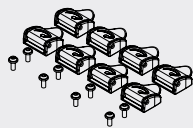9996098273

## General

We appreciate your choice of a new device from R82 - a worldwide supplier of technical aids and appliances for children and teenagers with special needs. To fully benefit from the options offered by this device please read this user manual before use and save it for future reference.

### Intended purpose

The family of R82 Standing devices consist of assistive devices intended for the alleviation of or compensation for a functional impairment due to an injury or disability. The devices are designed for an individual lacking the ability to maintain the standing position without the use of an assistive device.

Addax is an assistive device designed to support infants, children and adolescents (and adults who fit within the dimensions of the device) to perform supported standing. The device is intended for users where limited leg and trunk control makes standing upright without support difficult or impossible. The device has adjustable features in the frame and accessories to accommodate for the individual needs of the user.

### Intended user group

Selection of an appropriate device size is dependent on body weight and physical dimensions of the user. Secondary users of the device are caregivers providing assistance to the primary user as well as clinicians/technicians who set up the device.

The Addax is intended for infants, children and adolescents who require postural support to perform supported standing and activities with the upper extremities. It is designed for users with persistent, significant, or substantial functional and/or structural disabilities that affect their ability to stand independently. It is also suitable for adults who fit within the dimensions of the standing frame.

### Intended environment

The devices are intended for indoor use, in home environments or institutions. They are not to be used outdoor and in wet areas.

### Intended application

The user is the person standing in the device. The caregiver is the person operating the device.

The cumulative sum of time of use (exposure) will be either prolonged (more than 24 hours but less than 30 days) or long-term (exceeding 30 days). The use of the device will be discontinued multiple times during the day when the user will be lying down, sitting in another assistive device or a conventional product.

The device is intended for use with intact skin only. Depending on the configuration and the accessories used, the contact will often be with the skin of the head and neck region, shoulder region, upper extremities and lower extremities of the user and caregiver.

The device is intended for refurbishment and reuse.

### Indications

The assistive devices in the R82 Standing Family are indicated for children and adults who need postural support to enable supported standing due to a persistent, significant, or substantial impairment of mobility, i.e., disability. The indications for using the devices in the R82 Standing Family are based on functional ability and not a specific health condition or age.

The Addax is indicated for users with impairments in standing ability due to functional and/or structural damage to the legs or the legs and trunk muscles. This includes, but is not limited to, conditions such as neurological or neuromuscular diseases that result in postural instability and incorrect posture.

### Contraindications

There are no known contraindications.

### Precautions

Precautions for use are based on known and possible risks as identified in the R82 Standing devices User Manuals and detailed in the R82 Standing device Risk Management as follows.


#### Cautions:

Special attention should be given to the considerations about the configuration of the device and the use of accessories for users with limitation in their:

- Sensory functions and the sensation of pain, consider extra checks of the pressure caused by the support provided by the device and its accessories.
- Movement function (e.g. experience involuntary movements).
- Ability to maintain a standing position, use the device with accessories (e.g. chest or back support).
- Ability to maintain the position of their head, consider using accessories or to use a device that provide support for positioning of the head.

### Enclosed tools:

- 5 mm Allen key.
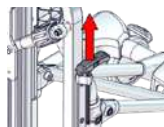
### CE Declaration of Conformity

This device conforms to the requirements of the Medical Devices Regulation (2017/745).

The CE mark must be removed when:

- Rebuilding the device
- Using the device outside the intended purpose and declaration of conformity
- Using other than original R82 spare parts and fittings.

### **Combinations**

When combining the device with other devices, the combination must not change the intended purpose of the devices or modify the devices in such a way that compliance with the applicable requirements may be affected. The person / company responsible for the combination of the devices must ensure that the applicable requirements are fulfilled.

Applicable requirements examples:

- Material safety (biocompatibility, flammability), mechanical safety (stability, squeezing hazards, inter connection), electrical safety.

### **Disposal**

When a device reaches the end of its lifetime, it should be separated by material type so that the parts can be recycled or disposed of properly. If required, please contact your local dealer for a precise description of each material. Please contact your local authorities to be informed of the options available for recycling in an environmentally friendly way.

### **Accessories and spare parts**

The devices from R82 can be supplied with a variety of accessories which comply with the needs of the individual user. Spare parts can be ordered on request. Find the specific accessories on our website or contact your local dealer for further information.

### **R82 Warranty**

[www.etac.com/support](http://www.etac.com/support)

### **Presales information and instructions**

Presales information and instructions are available at [www.etac.com](http://www.etac.com), under your device, in the documents section.

## Safety

Signs, symbols and instructions placed on this device must never be covered or removed and must remain present and clearly legible throughout the entire lifetime of the device. Immediately replace or repair illegible or damaged signs, symbols and instructions. Contact your local dealer for instructions.

Any adverse event in relation to the device must be reported to your local dealer and the national competent authority in a timely manner. The local dealer will forward information to manufacturer.

### Caregiver

- Read the user manual carefully before use and keep it for future reference. Incorrect use of this device may cause serious injury to the user.
- Settings and adjustments to the device or accessories are only to be made by persons who have read the user guide and who have been given the necessary instructions by a device advisor.
- This device contains small parts which could present a choking hazard if removed from their intended placing.
- Ensure that moving parts do not squeeze or trap body parts when adjusting them
- Ensure that body parts are not trapped in openings for mounting of supports/accessories.
- Always use the correct lifting techniques and aids for this purpose.
- Never leave the user unattended in this device. Ensure permanent supervision by a caregiver.
- Always use the necessary supports and fastenings to ensure that the user is safe and stable in the device.
- Repair/replacements must only be made using new original R82 spare parts and fittings and carried out in accordance with the guidelines and service intervals prescribed by the supplier.
- If there is any doubt as to the continued safe use of your R82 device or if any parts should fail, stop using the device immediately and contact your local dealer as soon as possible.
- Transportation of the device whilst the user is in it is not allowed. It is allowed to use the device with the user safely strapped in it for transfers over short distances. Do not cross any obstacles during transfer.
- Ensure that the user is supported, when performing any adjustment on the device or accessories.

### Environment

- Observe the surface temperature on the device, before placing a user in it. This goes especially for users with insensitive skin as they cannot feel the heat. If the surface temperature is above 41°C then allow the device to cool before use.
- Use this device on a level surface without any obstacles.

- When using on inclines, please refer to measurements.

### User

- If the user is close to the maximum user mass and/or has a lot of involuntary movements, e.g. rocking, a larger size with a higher maximum user mass or a different R82 device should be considered.

### Device

- Carry out all positional adjustments on the device and accessories and ensure that all knobs, screws and buckles are securely fastened before use. Keep all tools out of the reach of children.
- Ensure that all adjustable parts are securely fastened before use.
- Activate the brakes before placing the user into the device.
- Ensure that the stability of the device is maintained before placing the user in it.
- Do not expose the gas spring to pressure or high temperatures. DO NOT puncture.
- Do not expose the piston rod to scratches etc. It is important to use a cloth or the like, if tools are used to adjust the gas spring.
- This device is tested for resistance to ignition at level non-ignition in accordance with applicable standards.
- The device might lose its flame resistant characteristics if using other than original R82 cushioning.
- This device is not to be occupied in a motor vehicle. The user must be positioned in a vehicle seat and the device must be stored in a luggage compartment.
- The user must not occupy the device while it is being transported into or out of a vehicle by a lift.
- Unoccupied device is suitable for land and air transport.

## Service & maintenance information

Service and maintenance of a medical device is wholly the responsibility of the owner of that device. Failure to service and maintain a device in accordance with the instructions may invalidate the device's warranty. Furthermore, failure to service and maintain a device may compromise the clinical condition or safety of users and/or their caregivers. Do not perform service and maintenance while the user is in the device. Contact the local dealer for assistance, if needed, in setting up, using or how to perform service and maintenance of the device.

### Lifetime

The lifetime of this device in normal use is 8 years if all maintenance and servicing is carried out in accordance with the manufacturer's instructions and demonstrably recorded.

### Reuse

This device and its accessories/components are suitable for reuse.

### Refurbishment for reuse

Before reuse or changing the user, follow the service checklist, maintenance and washing instructions and carry out regular inspections.

### Service interval

The frequency of inspection may be altered depending upon use and wear. It is recommended that the device is inspected yearly, every time the device is re-issued for use and after a long-term storage (more than 4 months). The inspection should be carried out by a person who understands the use of the device.

### Service checklist:

Check and adjust the following as per their intended use:

- Cushioning.
- Nuts and bolts on the device are tightened.
- Fittings and buckles for signs of worn out parts.
- Inspect the device and all of its accessories and replace any worn out parts before use.
- Ensure the castors are fully-functional and secure.
- Ensure the foot supports are securely fastened.
- All labels on the device are intact.

Repair or replace any damaged or worn out parts.

### Maintenance

Before every use it is recommended to: Wipe off any residue and dirt from the device, using a cloth with warm water and a mild detergent/soap without chlorine and let it dry before use.

Carry out a visual inspection, for any damaged or worn-out parts.

### Washing

#### Main device

This device can be washed at 60 °C using a mild detergent for 10 minutes in a washing machine designed for medical devices. Use the machine's

drying function to dry the device. Unmount all accessories and wash them separately.

Tilt/rotate the device after washing to facilitate water drainage from within the device.

#### Fabric and covers

The cover can be washed at max. 60 °C using a mild detergent in a washing machine. Separate foam and cover before wash.

#### Foam

It is not recommended to wash the foam parts.

#### Handwash

This device and its accessories can also be washed by hand. Use warm water and a mild detergent/ soap without chlorine and let all parts dry before use.

#### Disinfection

The device can be disinfected with a 70% disinfectant IPA solution. It is recommended to wipe off any residue and dirt from the device, using a cloth with warm water and a mild detergent/soap without chlorine and let it dry before disinfection.

#### Materials

- Stainless steel
- Aluminium
- Plastic
- Rubber
- Polyurethane
- Neoprene
- Polyester
- Steel
- Brass

#### Surface treatment

The following surface treatment processes were used to protect against corrosion:

- Lacquered surfaces = Polyester powder coating or ED-coating
- Non-lacquered aluminium parts = Anodized coating
- Non-lacquered steel surfaces = Galvanized.

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**Etac Ltd.**
**Gelbart AG**