tauron|dts[®]

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Thank you for purchasing the tauron|dts shower commode chair from DIETZ.

This product is a medical aid. In order to use it properly and avoid possible risks due to improper use, you must familiarise yourself thoroughly with its operation.

Read this user manual carefully before you begin using your product. It contains important notes and information that will enable you to enjoy the full technical benefits. In addition, it will provide you with information to help keep your product in safe working order while ensuring it retains its value for as long as possible.

If you have any questions or require additional information, please contact the specialist retailer which supplied the product to you.

You can always find the latest information about your product on our website. For information or queries concerning product safety and on recalls, please contact DIETZ in writing or by telephone. You will find our contact information on the back of this user manual.

NOTE

A large-print version of the user manual is available for visually impaired persons at **shop.dietz-group.de** in PDF format.

1 Important information

1.1 General information	36
1.2 Intended purpose	36
1.3 Indications	36
1.4 Contraindications	37
1.5 Symbols	37

2 Product description

2.1 Scope of delivery	38
2.2 Accessories	38
2.3 Structure of the tauro dts standard side element	39
2.4 Structure of the tauro dts abduction side element with standing aid	39
2.5 Structure of the tauro dts featuring swivelling side element	39

3 Preparing the product for use

3.1 Assembly instructions	40
3.2 Tools	40
3.3 Seat width	40
3.4 Backrest height	42
3.5 Seat depth	42
3.6 Seat height	44

4 Using the product

4.1 tauro dts standard side element.	45
4.2 tauro dts abduction side element with standing aid	45
4.3 tauro dts swivelling side element.	45
4.4 Getting in and out.	46
4.5 Sideways transfers.	46
4.6 Using the product as a commode chair . . .	47
4.7 Using the product as a shower chair . . .	48
4.8 Rotating the seat surface.	48
4.9 Adjusting the backrest cushion	49
4.10 Folding the product.	49
4.11 Transport.	50
4.12 Accessories and add-ons from external suppliers	50

5 Technical specifications

5.1 DIETZ's measurement system.	51
5.2 Dimensions and technical information . .	52
5.3 Product tests.	52
5.4 Further information	53
5.5 Materials	53
5.6 Product labelling.	54

6 Safety instructions

6.1 Safety instructions.	55
----------------------------------	----

7 Care and maintenance

7.1 Servicing/maintenance.	56
7.2 Maintenance schedule	57
7.3 Troubleshooting.	58
7.4 Cleaning	59
7.5 Disinfection.	59
7.6 Disposal	59
7.7 Passing on the product for reuse	60
7.8 Storage.	60

8 Manufacturer's declarations

8.1 Warranty	61
8.2 Service life	61
8.3 Liability.	61

1.1 GENERAL INFORMATION

This user manual is included in the scope of delivery and forms an integral part of the product. The user should keep this manual to hand at all times and it must remain with the product if passed on to a third party for reuse. To check

that you are in possession of the latest version of the user manual, please contact DIETZ GmbH. Read the user manual before assembling and using the product.

1.2 INTENDED PURPOSE

The tauron|dts shower commode chair (max. user weight 250 kg) has been specially designed for people who are obese and have trouble showering and going to the toilet due to restricted mobility. The product is intended exclusively for use indoors on a firm and level surface.

The product is suitable for use in old-age and nursing homes, hospitals, rehabilitation centres,

day care centres and in private households.

It is not intended to serve as permanent seating and must not be used for this purpose.

If you have a custom-made device it will have been labelled as such. In this case, the intended purpose may differ from the one specified in this manual and you should refer to the documentation supplied with the product.

1.3 INDICATIONS

The use of this product is indicated in cases where persons with a physical disability cannot use or cannot easily access the existing shower or toilet facilities due to:

- ▼ Obesity
- ▼ Paralysis
- ▼ Loss of limb(s)/limb defects/limb deformities
- ▼ Joint contracture/joint damage
- ▼ Various illnesses/diseases

1.4 CONTRAINDICATIONS

The product must not be used in the case of:

- ▼ Perception disorders
- ▼ Severe balance disorders
- ▼ Inability to sit
- ▼ A user weight exceeding 250 kg
- ▼ Nor must it be used by children or anyone who is less than 146 cm tall and/or who weighs less than 40 kg (BMI < 17)

1.5 SYMBOLS

These symbols draw attention to passages of text that contain helpful advice for using and operating the product in everyday life.

WARNING

It is mandatory to observe and comply with warnings. They draw your attention to the fact that failure to comply may result in injury and/or damage to the product or the environment. _____

NOTE

Tips and advice which make it easier to use the various functions. _____

The product is partially pre-assembled at the factory and delivered inside a cardboard box. If possible, please keep the packaging in case you need to store the product inside it at a later date.

If the contents are noticeably damaged, inform the shipper immediately.

WARNING

Risk of poisoning: When burnt, packaging films can release toxic gases.

Risk of suffocation: Keep plastic bags out of the reach of children.

2.1 SCOPE OF DELIVERY

Upon receipt of the goods, please check immediately that the contents are complete and undamaged. The contents consist of:

- ▼ 2 side elements incl. 4 telescopic legs and 4 rubber feet
- ▼ 1 backrest incl. backrest cushion
- ▼ 1 seat cushion incl. 2 plug-in tubes, 2 toilet pan retaining rails, 4 mounting clips
- ▼ 1 tool set: 1 double open-end wrench, 1 allen key
- ▼ 1 bolt set: 2 size M6 x 50 mm hexagon socket head bolts, 2 knurled nuts
- ▼ 1 user manual

NOTE

The equipment features as well as the scope of delivery may vary depending on the country and market.

2.2 ACCESSORIES

- ▼ Round toilet pan incl. lid
- ▼ Rectangular toilet pan incl. lid

2.3 STRUCTURE OF THE tauro|n|dts STANDARD SIDE ELEMENT

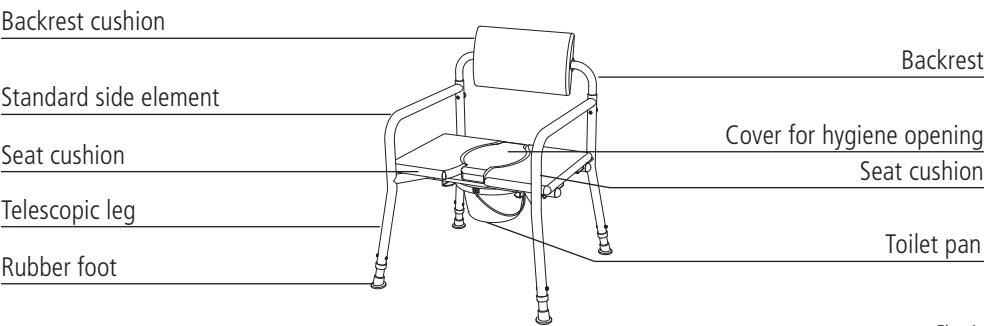
Fig. 1

2.4 STRUCTURE OF THE tauro|n|dts ABDUCTION SIDE ELEMENT WITH STANDING AID

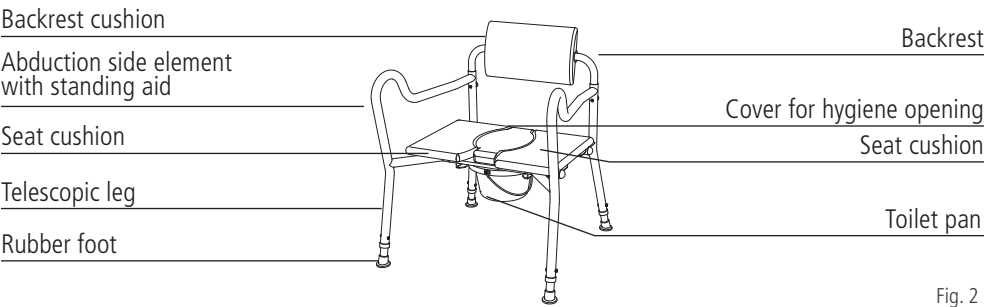
Fig. 2

2.5 STRUCTURE OF THE tauro|n|dts FEATURING SWIVELLING SIDE ELEMENT

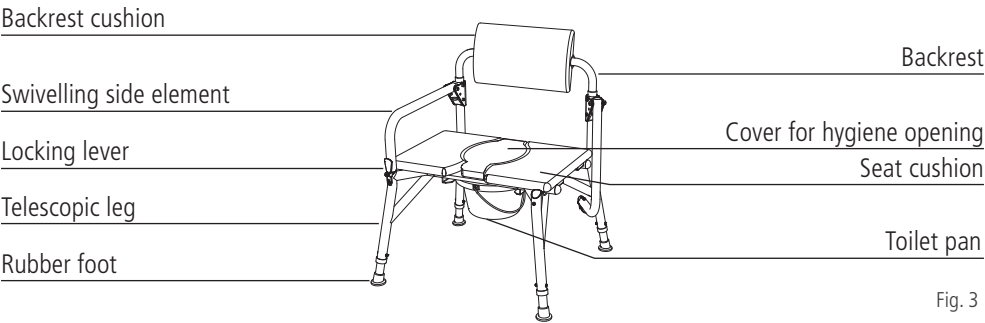
Fig. 3

3.1 ASSEMBLY INSTRUCTIONS

The assembly process and functional check are to be performed by the specialist retailer. The specialist retailer's technician will adjust the product as appropriate for the respective user.

NOTE

Initial set-up should be carried out by trained personnel of the specialist retailer

3.2 TOOLS

▼ Allen key (5 mm)

▼ Double open-end wrench (8 x 10 mm)

3.3 SEAT WIDTH

1. Remove the parts and bolt set from the cardboard packaging and place them all on a level, clean, sturdy surface (the floor or a table).
2. The seat can be set to one of three widths (560/630/700 mm). The backrest is pre-assembled for a seat width of 560 mm.
3. The first step is to slot the backrest into the side elements (1, Fig. 4). During assembly, take care to ensure that the parts are all aligned correctly and facing the right way. The side element sockets for the plug-in tubes should be pointing inwards. When slotting the backrest into the side elements, make sure that the backrest cushion is facing forwards.

WARNING

Risk of trapping: Be careful of your hands/fingers during assembly to avoid getting them trapped in moving parts.

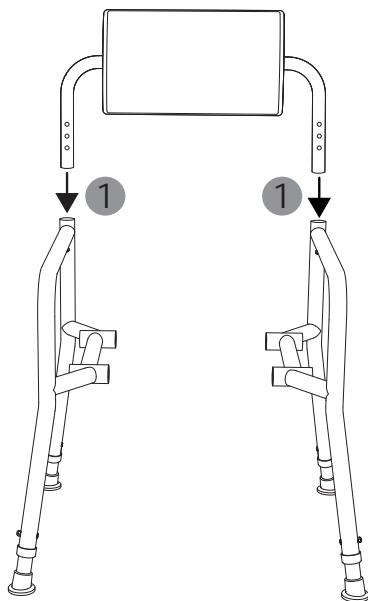
Fig. 4

4. Align the side elements so that they are parallel with each other (Fig. 5). Working from the front, insert a bolt (M6 x 50 mm) into each side element on the left and right (1, Fig. 5). Then secure them firmly with the knurled nuts on the rear of the product (2, Fig. 5).

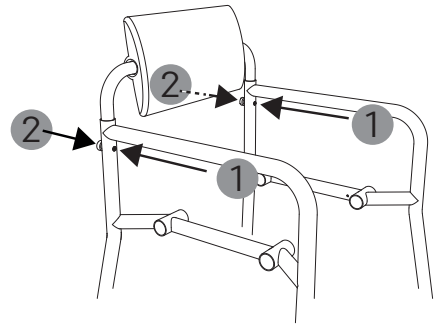
Fig. 5

5. To adjust the seat width so that it is right for you, proceed as follows: Loosen the knurled nuts on the rear of the backrest cushion (1, Fig. 6). Tilt the backrest cushion down towards the floor (2, Fig. 6).

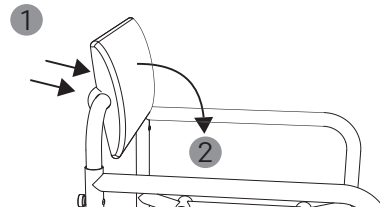
Fig. 6

6. The bolts on the rear of the cushion are now readily accessible and can be easily pulled out of the cushion (1, Fig. 7).

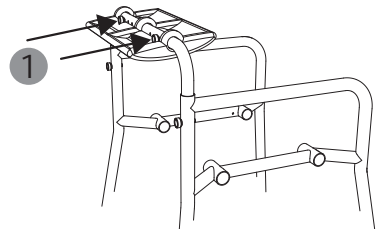
Fig. 7

7. The middle section inside the backrest is now free to move. There are three holes on each side for setting the required seat width (560/630/700 mm). To adjust the width, pull the side elements outwards (1, Fig. 8). When adjusting the seat width, make sure that the backrest remains centrally aligned.

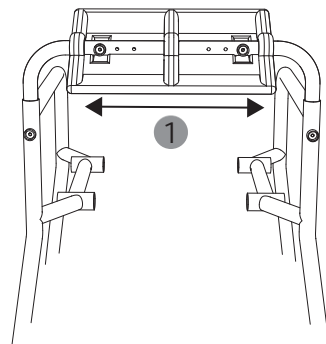
Fig. 8

8. Once you have selected the required seat width, reinsert the bolts and secure them with the knurled nuts.

3.4 BACKREST HEIGHT

1. To adjust the backrest height, loosen the knurled nuts on the side elements and pull out the bolts (1, Fig. 9).
2. There are three holes for adjusting the height of the backrest (410/435/460 mm). To adjust the height, pull both sides of the backrest up evenly (2, Fig. 9).
3. Once you have set it to the required height, secure the backrest again with the bolts and knurled nuts. Make sure you insert the bolts from the front (1, Fig. 10) and that you fasten the knurled nuts securely on the rear of the product (2, Fig. 10).

WARNING

Risk of trapping: Be careful of your hands/fingers during assembly.

Risk of accident: The product is only to be used with the backrest securely attached. _____

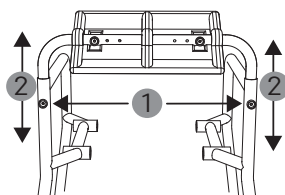
Fig. 9

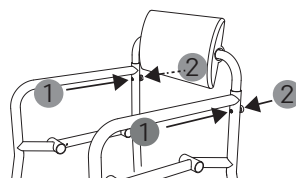
Fig. 10

3.5 SEAT DEPTH

You can mount the seat surface in one of two positions (450/510 mm). The seat depth is measured from the backrest tube to the front edge of the seat surface (see Chapter 5.1).

1. Select the seat depth that is right for you (Fig. 11, 12). Please note that the selected seat depth also determines where the hygiene opening is located with the cover attached to the seat surface.


WARNING

Risk of trapping: Be careful of your hands/fingers during assembly to avoid getting them trapped in moving parts.

Risk of accident: If you are able to lift the seat surface up at either the front or the back following assembly, it means that the plug-in tubes have not been assembled correctly.


Risk of accident: The plug-in tubes must be fully inserted and locked with the mounting clips before use. If the mounting clips are inserted from the front/rear, they may come loose. _____

2. Place the seat surface centrally on the tube sockets between the side elements, ensuring that it is lying straight (1, Fig. 13). To check that it is fitted correctly, shake the seat surface gently. If the seat surface no longer moves, it means that the locating bracket underneath the seat surface has been positioned correctly.
3. Working from the side (left or right), insert both plug-in tubes into the sockets underneath the seat surface (if necessary, gently twist the plug-in tubes to make them slide in more easily). Please ensure that the plug-in tubes are fully inserted (2, Fig. 13).
4. Try lifting up the seat surface to check that it has been securely mounted. You should not be able to lift it off the tube sockets at either the front or the back.
5. Twist the plug-in tubes so that the openings of the holes are pointing upwards and downwards. Insert the four mounting clips from below (1, Fig. 14).
6. The mounting clips prevent the plug-in tubes from slipping out sideways.

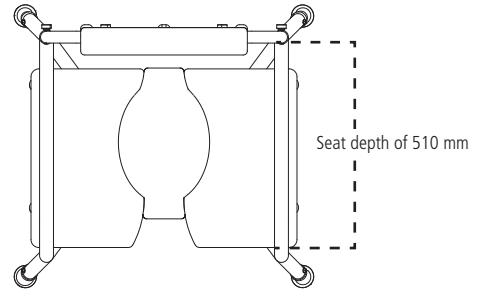Seat depth of 510 mm, hygiene opening at the front

Fig. 11

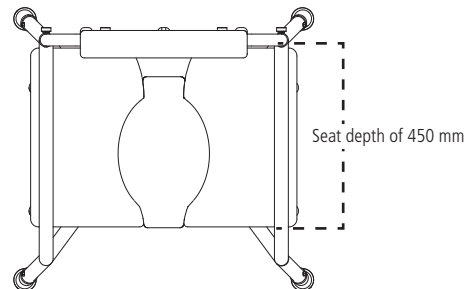Seat depth of 450 mm, hygiene opening at the rear

Fig. 12

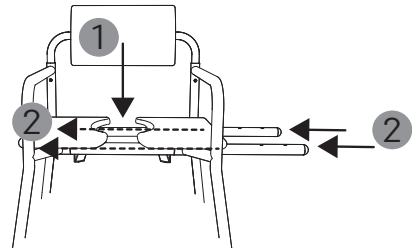
Fig. 13

NOTE

If the plug-in tubes cannot be fully inserted, proceed as follows: Loosen all four knurled nuts on the backrest and gently shake the side elements. You should now find that the plug-in tubes can be inserted with ease. Finally, retighten all the knurled nuts.

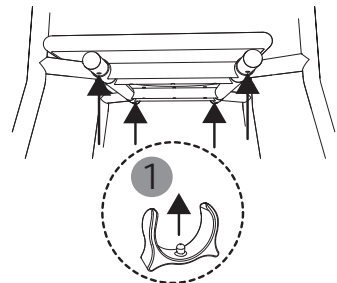
Fig. 14

3.6 SEAT HEIGHT

1. Use the tool set included in the scope of delivery to adjust the seat height.
2. To adjust the telescopic legs, undo and remove the bolts (1, Fig. 15).
3. You can now set the telescopic legs to the required length (2, Fig. 15), ensuring that the holes of the inner and outer tubes are aligned with each other.
4. Finally, reinsert the bolts, attach the nuts externally and tighten them.

WARNING

Risk of trapping: Be careful of your fingers when adjusting the telescopic legs.

Risk of tipping: Set all 4 telescopic legs to the same length.

Risk of accident: Make sure the bolts are firmly seated.

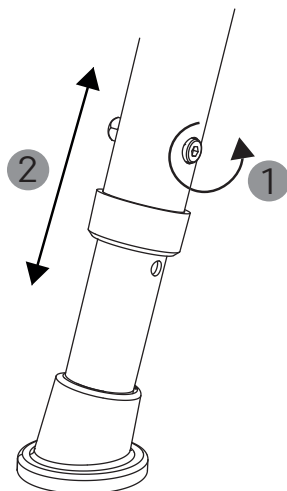
Fig. 15

4.1 TAURON|DTS STANDARD SIDE ELEMENT

These side elements (Fig. 16) are designed to meet the general requirements of everyday care settings.

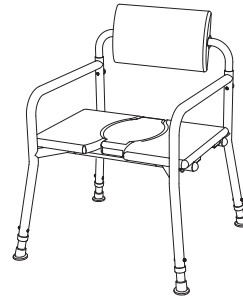
Fig. 16

4.2 TAURON|DTS ABDUCTION SIDE ELEMENT WITH STANDING AID

The abduction side elements (Fig. 17) add an extra 150 mm onto the seat width so that those who need considerably more space towards the front can sit comfortably on the chair.

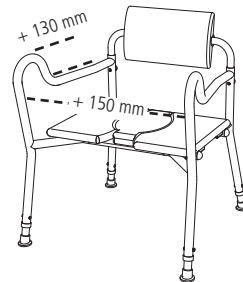
Fig. 17

4.3 TAURON|DTS SWIVELLING SIDE ELEMENT

The swivelling side elements on both sides (Fig. 18) enable easy transfer, e.g. from a wheelchair into the shower commode chair, and vice versa.

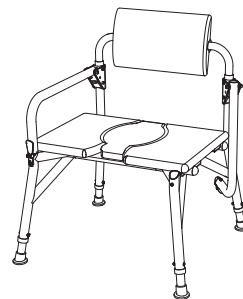
Fig. 18

4.4 GETTING IN AND OUT

To get in and out of the chair as safely and easily as possible, please proceed as follows.

Getting into your chair:

1. Position the product with all four legs standing securely on the floor and make sure there is nothing to hinder its use.
2. Position yourself in front of the chair, standing as close as possible to it. Turn around (Fig. 19) and carefully move backwards until you feel your legs touch the edge of the seat surface.
3. Grasp the side elements with both hands and use them to support yourself
4. as you safely and slowly lower yourself into the shower commode chair.

Getting out of your chair:

To get back out of the shower commode chair, follow the above steps in reverse.

WARNING

Risk of trapping! When inserting the cover for the hygiene opening, make sure it is inserted so that it lies straight.

Risk of tipping: Please note that body movements shift the centre of gravity and can increase the risk of the product tipping over.

Risk of tipping: Support yourself vertically on both armrests and not in any other way.

Risk of accident: Do not stand with your feet on the seat surface.

Risk of slipping: Do not allow any oily or highly soapy substances to come into contact with the rubber feet.

Fig. 19

4.5 SIDEWAYS TRANSFERS

Sideways transfers are only possible if you have opted for the swivelling side element model.

1. Position the shower commode chair at the point where you want to transfer yourself into or out of it.
2. To swivel the side element, release the locking lever by pulling it towards you (1, Fig. 20).

WARNING

Risk of accident: The armrest must always be swivelled down for transfers. Following transfer, lock the armrest again. Never lean on the side element or use it to support yourself while it is unlocked.

Risk of trapping: Be careful of your legs, fingers, hands and – where applicable – your clothing when swivelling the side element.

3. Release the side element from its fastener by lifting it up and to the side (2, Fig. 20).
4. The side element can now be swivelled down (3, Fig. 20).
5. You can now transfer yourself into or out of the chair by moving sideways.
6. After the sideways transfer, lock the side element again by following the above steps in reverse. The side element audibly clicks into place.

We recommend that you get an assistant to help you when transferring yourself out of the chair.

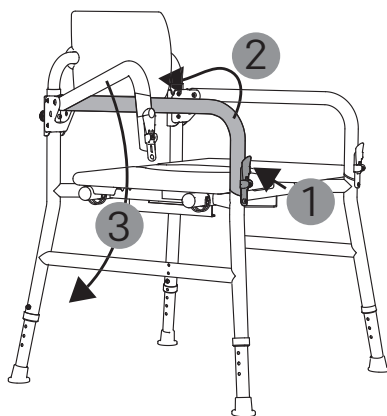
Fig. 20

4.6 USING THE PRODUCT AS A COMMODE CHAIR

Before use:

1. Without putting the lid on, slide the toilet pan over the guide rail under the seat from the rear/front so that it is positioned between the plug-in tubes (1, Fig. 21).
2. Remove the cover from the hygiene opening (2, Fig. 21).

WARNING

Risk of trapping: Do not insert the pan or pull it out while a person is still sitting on the chair.

Risk of infection: The pan must be emptied and cleaned after use. _____

After use:

1. To empty the pan, carefully pull it out towards the rear/front (1, Fig. 21).
2. After emptying it, clean/disinfect the pan and seat surface in accordance with the instructions.
3. Then reinsert the cover for the hygiene opening, ensuring that it is lying straight.

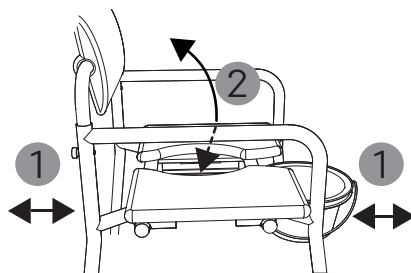
Fig. 21

4.7 USING THE PRODUCT AS A SHOWER CHAIR

Before use:

1. The toilet pan should not be inserted when using the product to take a shower. If necessary, carefully remove it by pulling it out towards the rear/front (1, Fig. 21).
2. It is entirely up to you whether you insert the cover for the hygiene opening into the seat surface or remove it (2, Fig. 21).

After use:

1. Always clean/disinfect the product in accordance with the instructions.

WARNING

Risk of trapping! Do not insert the pan or pull it out while a person is still sitting on the chair.

Risk of trapping! When inserting the cover for the hygiene opening, make sure it is inserted so that it lies straight.

4.8 ROTATING THE SEAT SURFACE

1. Remove the four mounting clips from the plug-in tubes (1, Fig. 22).
2. Next, push the plug-in tubes out sideways (left or right) by pressing the end caps (1, Fig. 23) so that the tubes protrude on the opposite side. You can now pull out the plug-in tubes (1, Fig. 23.)

WARNING

Risk of trapping: Be careful of your hands/fingers when assembling the seat surface, and when pushing in and pulling out the plug-in tubes.

Risk of accident: The mounting clips must be inserted before use. Make sure that the mounting clips enclose the plug-in tubes tightly, otherwise they may come loose.

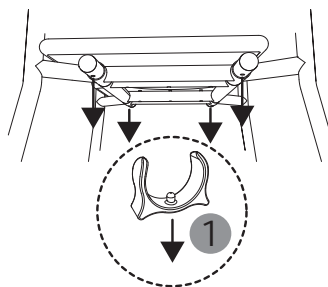
Fig. 22

3. The seat surface is now free to move and can be raised (2, Fig. 23) and rotated 180° (3, Fig. 23).
4. To secure the seat surface, please follow the instructions in Chapter 3.5, starting from no. 2.

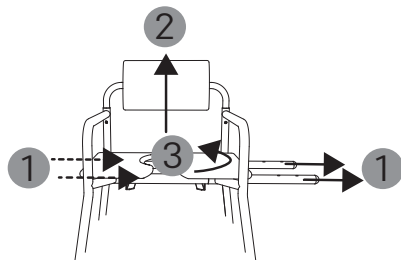
Fig. 23

4.9 ADJUSTING THE BACKREST CUSHION

The backrest cushion can be rotated 360°. It automatically adapts to your individual sitting position (Fig. 24).

NOTE

Positioning: Always align the backrest cushion to match how you are seated. _____

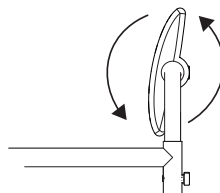
Fig. 24

4.10 FOLDING THE PRODUCT

1. Release the seat surface (Chapter 4.8, nos. 1. and 2.) and remove it from the product.
2. Undo the knurled nuts on the rear of the backrest and remove the bolts (1, 2, Fig. 25). Keep the bolts somewhere safe until you are ready to reassemble the product.

WARNING

Risk of trapping: Be careful of your hands/fingers when folding the product. _____

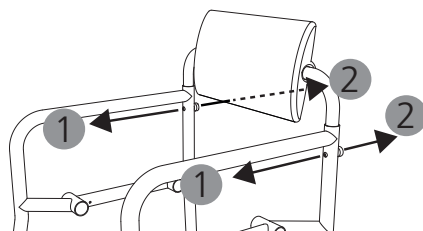
Fig. 25

3. You can now fold the product easily (1, 2, Fig. 26) by swivelling one of the side elements inwards and the other one outwards.
4. If you want to start using the product again, unfold it by performing the above steps in reverse and then follow the instructions in Chapter 3.3, starting from no. 4.

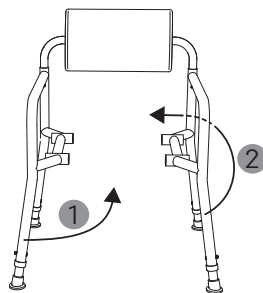
Fig. 26

4.11 TRANSPORT

The shower commode chair is not designed for frequent transport. Should it be necessary to transport the chair, we recommend folding it first (Chapter 4.10). The chair should be packed so that it is standing in a stable position and there is no danger of the means of transport performing sudden and dangerous movements. If you need to carry the product, we recommend lifting the chair by both side elements.

4.12 ACCESSORIES AND ADD-ONS FROM EXTERNAL SUPPLIERS

As a general principle, only original accessories from DIETZ GmbH are to be used. If third-party products are installed on the product, responsibility for the safety of the product passes to the person who installs the accessories or fits the add-on. In this case, the installer or fitter must issue a new declaration of conformity for the combination of accessories and product or add-on and product. The conformity declared by DIETZ GmbH according to MDR 2017/745, Annex II is rendered null and void.

WARNING

Safety hazard: The safety of the product cannot be guaranteed if add-ons/accessories are sourced from anyone other than DIETZ.

If accessories or add-ons are fitted to the product, then the safety instructions in the user manual for the accessories or add-ons must always be observed. _____

5.1 DIETZ'S MEASUREMENT SYSTEM

The dimensions given here refer to the standard configuration of the product.

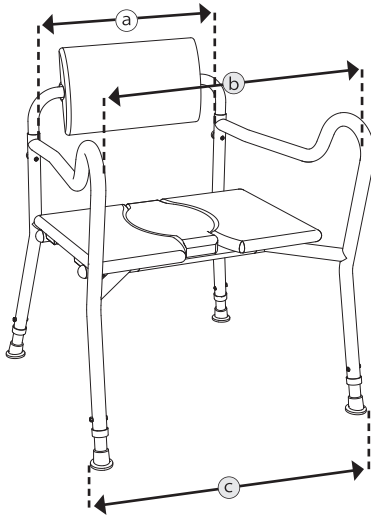
Fig. 27

NOTE

All dimensions measured on textile parts (e.g. backrest height) are to be understood as approximate values that are subject to a tolerance of ± 10 mm.

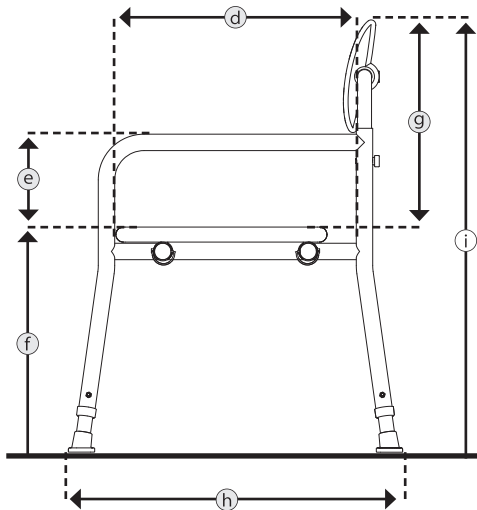
Fig. 28

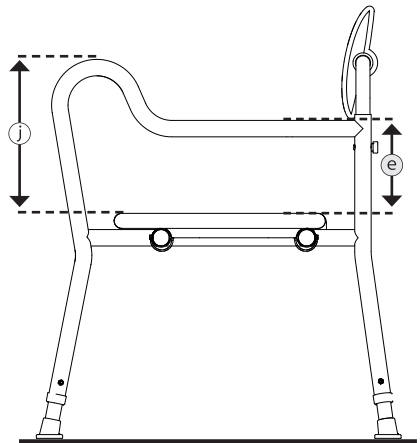
Fig. 29

5.2 DIMENSIONS AND TECHNICAL INFORMATION

Specifications			Standard	Abduction with standing aid	Swivelling
a	Seat width ¹	mm	560 / 630 / 700		
b	Seat width ¹ at front	mm	–	710 / 780 / 850	–
d	Seat depth ²	mm	450 / 510		
f	Seat height	mm	450 – 600		
h	Total length ³	mm	705 – 747	770 – 806	705 – 747
c	Total width ³	mm	SB 560: 750 – 792 SB 630: 820 – 862 SB 700: 890 – 932	SB 560: 810 SB 630: 885 SB 700: 960	SB 560: 750 – 792 SB 630: 820 – 862 SB 700: 890 – 932
i	Total height	mm	860 – 1060 (in 25 mm increments)		
	Seat surface (W × D × H)	mm	670 × 450 × 32		
	Backrest cushion (W × D × H)	mm	420 × 70 × 245		
	Adjustable backrest angle	°	360		
g	Backrest height	mm	410 / 435 / 460		
e	Side element height	mm	200		
j	Side element height at front	mm	–	330	–
	Armrest length	mm	580	680	580
	Drive-over height	mm	343 – 493		
	Round toilet pan (W × H / diam.)	mm	300 × 160 / 300		
	Rectangular toilet pan (L × W × H)	mm	410 × 300 × 160		
	Mass of heaviest part ⁴	kg	6,25	6,75	7,25
	Total weight	kg	14	14.5	15
	Maximum user weight ⁵	kg	250		

5.3 PRODUCT TESTS

ISO 17966	Tested according to the requirements and test methods for assistive products for personal hygiene
EN 1021-2	Tested for resistance to ignition by a flame in the case of upholstery materials and plastic components
ISO 10993-5	Tested for toxicity in the material
Test dummy weight	250 kg

¹ Measured between the side elements.² Measured from the backrest tube to the front edge of the seat surface.³ Depending on the seat height.⁴ Weight of shower commode chair excluding removable parts (in this case: without seat surface, plug-in tubes and toilet pan).⁵ Maximum user weight including any additional load. The weight of the additional load reduces the maximum user weight accordingly.

5.4 FURTHER INFORMATION

Specifications		tauron dts
Colours		Mangan Grey-01, Mangan Grey-02, Lime Green
Standard side element	kg	1.4
Abduction side element with standing aid	kg	1.6
Swivelling side element	kg	1.9
Backrest incl. backrest cushion	kg	1.45
Seat cushion	kg	5.2
Plug-in tube	kg	1.2
Round toilet pan	kg	0.54
Rectangular toilet pan	kg	0.2
Telescopic leg	kg	0.17
Rubber feet	kg	0.05
Ambient temperature during operation		+10°C to +40°C.
Ambient conditions during storage		0°C to +45°C // 20% to 75% relative humidity

5.5 MATERIALS

Side element	Aluminium, powder coated
Backrest	Aluminium, powder coated
Backrest cushion	PU (polyurethane)
Knurled nuts	PA (polyamide), brass
Seat cushion	PU (polyurethane), wood
Locating bracket, seat	PP + 15% GF (polypropylene + glass fibres)
Plug-in tubes	Stainless steel
End caps for plug-in tubes	PP (polypropylene)
Mounting clip	PP (polypropylene)
Round toilet pan	PP (polypropylene)
Rectangular toilet pan	PP (polypropylene)
Toilet pan retaining rail	Aluminium, powder coated
Telescopic leg	Aluminium, powder coated, PP (polypropylene)
Rubber feet	TPE (thermoplastic elastomer), stainless steel

Note: All metals used are corrosion-resistant.

5.6 PRODUCT LABELLING

The product label and warning signs on the product must remain legible. Please have any illegible or missing labels/signs replaced by your specialist retailer immediately. The serial number is very important for identifying the product and can be found on the product label. It must not be removed.

■ Labelling for custom-made devices

SONDERANFERTIGUNG
CUSTOM MADE DEVICE

Fig. 30

Custom-made devices are products which have been made specifically for a particular user. In this case, the product label, warnings on the product and possibly also the intended purpose will differ compared to a standard device. Please pay attention to the additional documentation.

■ Product label

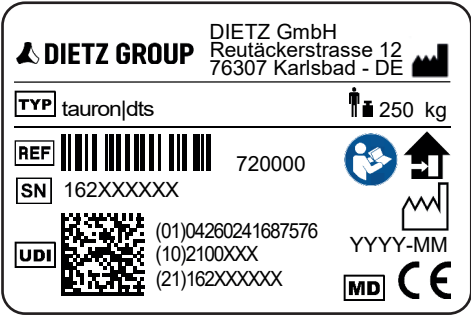
Fig. 31

TYP	TYP = Designation/model
REF	REF = Item number
SN	SN = Serial number
UDI	UDI = Unique Device Identifier (01) UDI-DI/GTIN (10) Batch/lot number (21) Serial number
	Manufacturer/manufacturer's logo/contact address
	Maximum user weight incl. additional load
	Date of manufacture
	The product is only approved for indoor use.
	Follow the user manual.
MD	Medical Device
CE	CE marking MDR 2017/745 Annex VIII

6.1 SAFETY INSTRUCTIONS

When using the product, always observe the following safety instructions to prevent falls, hazardous situations and damage to the product:

- ▼ The product is only to be used indoors and on a firm, level surface.
- ▼ Do not use the product as a climbing aid or stand on it.
- ▼ Do not load the product unevenly.
- ▼ There is a risk of slipping when getting in and out of the chair in rooms with wet floors, particularly after shower products (e.g. soap, liquid soap or shampoo) have been used.
- ▼ Check the product for damage before each use. Do not use the product if damaged.
- ▼ After disassembling and reassembling the product, always check its stability and ensure that all telescopic legs are set to the same length.
- ▼ If your mobility is severely restricted, get a care assistant to help you when using the product.
- ▼ The product must not be used in water tanks or swimming pools.
- ▼ The surface of the product may become hot due to hot water or exposure to sunlight. Risk of burns!
- ▼ Serious incidents that occur in relation to the product must be reported to the manufacturer and the responsible authorities.
- ▼ If your product is custom-made (see product labelling), the additional documentation – including all safety instructions – must be observed in addition to the latest version of the user manual.

7.1 SERVICING/MAINTENANCE

We recommend performing the inspections listed in the maintenance schedule at regular intervals to ensure that the product is always safe to use. Improper or insufficient product care and maintenance will result in restricted liability on the part of the manufacturer. The maintenance schedule does not provide any information on the amount of work that will actually need to be performed on the product.

NOTE

As the user, you will be the first to notice any damage. If you notice a fault as described in the maintenance schedule or further defects and functional impairments, contact an authorised specialist retailer immediately. _____

WARNING

To maintain the operational safety of the product, repairs must only be performed by specialist shops using DIETZ original spare parts. The relevant spare parts catalogue can be found on our website by looking under the information for the product concerned. _____

NOTE

Unless otherwise stated, the inspections and measures stipulated in the maintenance schedule must be performed by the user or their care assistant. _____

NOTE FOR SPECIALIST RETAILERS

If the user reports product abnormalities, please check all the inspection points listed in the maintenance schedule. The checks must also be performed after longer storage periods (>4 months) and whenever the product is going to be reused. _____

7.2 MAINTENANCE SCHEDULE

What	Description	Monthly	
		Before use	
Seat surface and backrest cushion Check for damage	- The seat surface and backrest cushion must not be damaged. - In the event of functional or surface damage, the components must be replaced by a specialist retailer.	x	
Adjustable parts Check for firm seating	Check that all screws/bolts for adjusting the product are firmly seated.	x	
Plug-in tubes Check for correct function and wear	Check that the mounting clips are seated correctly. They must be inserted from below and tightly enclose the plug-in tubes.	x	
	- Check the plug-in tubes to make sure they are not bent.  Tubes that are bent must be replaced by a specialist retailer.		x
Swivelling side element Check for correct function and wear, and for firm seating.	Check that the side elements are securely locked.	x	
	- Check the locking mechanism for functional damage.  If damage of this type is found, discontinue use of the product, as safe use can no longer be guaranteed.		x
Toilet pan Check for correct function and wear	- The pan should be easy to slide in and out. - In the event of functional or surface damage (toilet pan, guide rail), the components must be replaced by a specialist retailer.	x	
Bolted connections Check for firm seating	- All fastening bolts must be firmly tightened.  Repeated loosening and re-tightening causes self-locking nuts and bolts to lose their effectiveness. If this happens, they must be replaced by an authorised specialist retailer.		x
Frame, complete: Proper functioning and damage	- The frame parts (side elements, backrest, sockets for plug-in tubes) must not show any signs of damage, cracks or bending.  If damage of this type is found, discontinue use of the product, as safe use can no longer be guaranteed.		x
Visual inspection Check for loose parts, breaks, corrosion or other damage	 If damage of this type is found, discontinue use of the product, as safe use can no longer be guaranteed.		x
Entire product Check for dirt	The seat surface/backrest cushion and toilet pan must be cleaned and disinfected after every use.	x	
	The entire product must be cleaned based on how dirty it is, but at least once a month (see "Cleaning and disinfection").		x

7.3 TROUBLESHOOTING

What	Possible causes	Measures
The seat surface is not straight.	• The telescopic legs have been set to different lengths. • The plug-in tubes are not positioned correctly (Chapter 3.5 Seat depth).	Correct the setting.
The shower commode chair wobbles.	• The telescopic legs have been set to different lengths. • Check that all knurled nuts are firmly seated.	Correct the setting.
The cover for the hygiene opening cannot be inserted.	• The correct procedure for inserting the cover has not been followed.	Rotate the cover 180° and try inserting it that way round.
It is not possible to slide the toilet pan in/out.	• The guide rail has come loose on one side.	Tighten the bolts.
	• The toilet pan is damaged. • The guide rail is bent.	Consult an authorised specialist retailer.
It is not possible to lock/unlock the swivelling side element.	• The locking receptacle is not free or the locking lever is dirty.	Clean the components.
	• The locking lever is defective.	Consult an authorised specialist retailer.
Squeaking and grinding noises can be heard.	• There are various possible causes.	Consult an authorised specialist retailer.

7.4 CLEANING

To clean the product, use lukewarm water mixed with a mild detergent (without solvents and with a pH value of 5 to 9).

After cleaning, wipe with clear water and then dry the chair with a cloth.

WARNING

Material damage: Do not use any of the following for cleaning purposes: high-pressure cleaners; strong, corrosive chemicals; wire brushes; or scouring agents.

Risk of corrosion: To avoid corrosion, ensure that all materials are completely dry after cleaning.

7.5 DISINFECTION

Wipe down or spray all parts thoroughly with a commercially available surface disinfectant. Take particular care to disinfect surfaces that come into frequent contact with hands/skin, such as the seat surface, toilet pan, backrest incl. cushion, and side elements.

Please only use the following disinfectants; use of other chemicals might damage the product.

- ▼ Aldehyde-free disinfectants based on alcohol (max. 70% propyl alcohol)

WARNING

Material damage: During disinfection, observe the use and application instructions issued by the cleaning agent/disinfectant manufacturer.

Risk of infection: The product must be disinfected before each change of user. This is important for avoiding cross-contamination.

7.6 DISPOSAL

If you no longer need the product, please contact your specialist retailer. They will then pick up the product and either dispose of it properly or make further use of it. Otherwise, please take the product to your local waste disposal facility.

7.7 PASSING ON THE PRODUCT FOR REUSE

Your product is suitable for transfer to/reuse by others, unless it is a custom-made device and marked accordingly. The number of times it can be reused depends on the extent of material wear and whether the product is still in good working order.

When passing on the product for reuse, please remember to hand over all the technical documents necessary for safe handling and use to the new user or the specialist retailer.

Before it is reused, the product must be cleaned, disinfected, checked for damage and approved by a specialist retailer.

As part of this, all the inspection points listed in the maintenance schedule must be checked on the product.

7.8 STORAGE

If you wish to store the product, ensure that it is kept in a dry place protected against direct, strong sunlight and at a temperature of 0°C to +45°C. To save space during storage, you can fold the product (Chapter 4.10).

After longer periods of storage (> 4 months), please check all the inspection points on the product as listed in the maintenance schedule before reusing/redeploying the product.

WARNING

Risk of infection: The product must be disinfected before each change of user. This is important for avoiding cross-contamination. _____

NOTE

Products marked as custom-made are excluded from transfer and reuse. _____

WARNING

Do not store the product near a heat source and do not place any objects on the product during storage. _____

8.1 WARRANTY

The warranty covers all product defects which are verifiably attributable to material or manufacturing flaws. The warranty period is 24 months as of receipt of notification of readiness for shipping, but following delivery at the latest. The warranty does not cover damage which has

occurred due to natural wear and tear, or has been caused deliberately or through negligent or improper operation or use. This also applies to the use of unsuitable care products, lubricating oils or greases.

8.2 SERVICE LIFE

The product has an expected service life of five years when used daily and exactly as intended. This assumes compliance with the maintenance

and safety specifications stipulated in this manual. The indicated service life does not constitute any kind of additional warranty.

8.3 LIABILITY

DIETZ GmbH only bears liability if the products are used under the specified conditions and for the specified purposes. It is very important to handle and care for the products properly as per the instructions. DIETZ GmbH is not liable for damage caused by components and spare parts that have not been approved by DIETZ GmbH.

Repairs must only be carried out by authorised specialist retailers or by the manufacturer itself.

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© DIETZ GmbH, Karlsbad
Nachdruck, auch auszugsweise, nur mit schriftlicher Genehmigung von DIETZ GmbH, Karlsbad.
Not to be reproduced, whether in part or in whole, without the written permission of DIETZ GmbH, Karlsbad.
Verveelvoudiging, ook deels, uitsluitend met schriftelijke toestemming van DIETZ GmbH, Karlsbad.
La reproduction, même partielle, n'est autorisée qu'avec l'accord écrit de DIETZ GmbH, Karlsbad.
Riproduzione, anche parziale, solo previa autorizzazione scritta di DIETZ GmbH, Karlsbad.